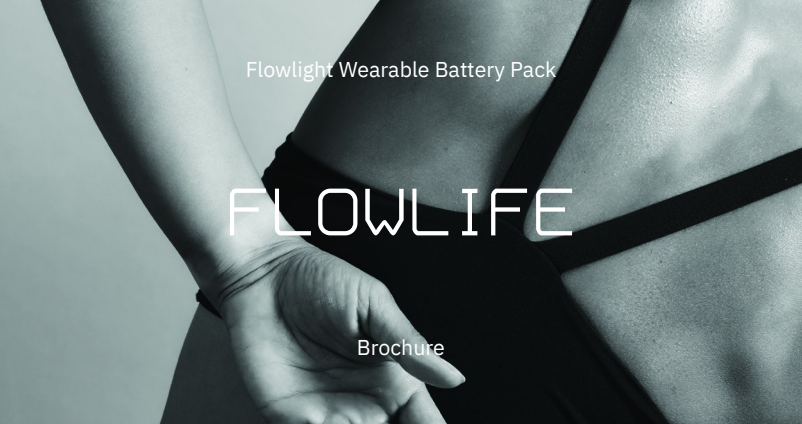
Flowlight Wearable Battery Pack

FLOWLIFE

Brochure

Flowlight Wearable Battery Pack

All-in-One Control Unit for the Flowlight Wearable Recovery System

Table of Content

4 — Manifest

6 — Flowlight Wearable Battery Pack

7 — Accessories

8 — Flowlight Wearable Battery pack: A Comprehensive Guide

14 — Starting your Journey

21 — Contact Information

24 — Manual

- | | | |
|-----------|------------|---------|
| • English | • Deutsch | • Dutch |
| • Svenska | • Français | • Suomi |
| • Norsk | • Espanol | |
| • Dansk | • Italiano | |

Manifest

In a world where limits are meant to be shattered, we fuel the fire of relentless pioneers. Our brand champions those who break the mold and embrace the extraordinary, symbolizing a life lived without compromise.

Committed to relentless innovation, we design for those who demand excellence and performance, understanding that the journey into the unknown often leads to unforgettable achievements.

We want to power dreams, aspirations, and triumphs of the audacious, sparking the passions of the unstoppable.


For those who redefine boundaries and explore new horizons, we are Flowlife. Made Possible.

Why Flowlight Wearable Battery pack?

The Flowlight Wearable Battery pack is the control centre and power source for your Flowlight wearable therapy device — all in one compact unit. It magnetically clips directly onto the wearable, giving you instant control over every session through a clear LCD screen and five dedicated buttons. Charge it once via USB-C, clip it on, and your recovery session is ready.

- **One unit. Full control** - The Flowlight Wearable Battery pack combines session control and built-in battery power into a single magnetic unit. Press to turn on, choose your mode and intensity, set your timer, and start — directly from the controller clipped to your wearable.
- **Magnetic clip. Instant attachment** - Magnetic connectors attach the controller firmly to the wearable — no fumbling with cables or separate power packs. The secure connection maintains full power delivery throughout your session.
- **LCD display. Five-button precision** - The backlit LCD screen shows your current mode, intensity level, and session timer in real time. Five intuitive buttons — Power, Mode, Start/Pause, +, and — give you full control without lifting the device or opening an app.

Accessories

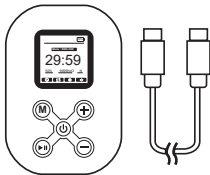What's Included

Flowlight Wearable Battery pack
(control unit with built-in battery)

EU power adapter

Type-C USB charging cable

User manual

Works With Flowlife Wearables

- Flowlight Wearable Knee 60 Two Waves
- Flowlight Wearable Foot 60 Two Waves

Flowlight Wearable Battery pack: A Comprehensive Guide

The Flowlight Wearable Battery pack attaches magnetically to your Flowlight wearable and provides both the power supply and full session control through its built-in LCD display and five-button interface. It charges via USB-C with an integrated charging circuit — no external adapter required.

Why it matters:

- Magnetic attachment keeps the controller in place throughout movement and therapy sessions
- LCD display provides real-time feedback on mode, intensity, and remaining session time
- Five dedicated buttons give precise, tactile control — no app required
- Built-in 6000mAh battery delivers multiple full sessions before recharging
- USB-C charging with built-in circuitry — simply connect to any USB-C power source

Magnetic Connection

Align the magnetic connectors on the Flowlight Wearable Battery pack with the corresponding connection points on your Flowlight Wearable. The magnets will seat firmly and automatically — you will feel a secure click. A secure connection is confirmed when the LCD screen activates on power-on.

LCD Display

The LCD screen displays current mode (Red / NIR / Combined / Pulse), intensity level, and session timer. During a session, the display updates in real time so you can monitor progress without removing the device.

Five-Button Control

POWER 

Hold to turn the controller on or off. Mode: cycles through settings — press once to select light mode, again to set intensity, again to set Pulse frequency. session timer from the home screen

+ and -  

Adjust the current setting up or down, or set the session timer from the home screen.

Start/Pause 

Begins the session or pauses it mid-session.

Mode

Cycles through settings. Press to step between light mode, intensity, and pulse frequency.


Charging the Battery Pack

Connect the included USB-C cable to the charging port on the controller. The LCD display or indicator light will show charging status. The LED indicator shows 100% when fully charged. Charging time is approximately 2.5 hours from empty.

Starting
Your Journey

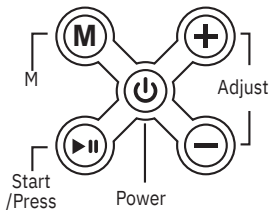1. Charge the Flowlight Wearable Battery pack fully before first use — connect via USB-C until the indicator shows 100%.
2. Fit your Flowlight Wearable onto the target area and secure the straps.
3. Align the Battery Pack magnetic connectors with the wearable connection points — the magnets will seat firmly.
4. Hold the Power button to turn the controller on — the LCD screen will illuminate.
5. Press Mode and use + / - to select your light mode: Red / NIR / Combined / Pulse.

6. Press Mode again and use + / – to adjust intensity level.
7. Press Mode again and use + / – to set your Hz level.
8. In the home screen view, press + / – to set your session timer.
9. Press Start to begin your session — the timer will count down on the LCD display.
10. Press Start/Pause at any time to pause mid-session. Hold Power to end the session early.

Targeted Joint Recovery

Precise therapy.
Hands-free
control.

Recovering joints need consistent, correctly dosed light therapy — but fumbling with settings mid-session or losing the controller breaks the routine and reduces effectiveness.

Magnetic
Control Unit

The Flowlight Wearable Battery pack clips directly onto the wearable and stays in place during use. The LCD screen and five buttons let you adjust mode and intensity in seconds — without detaching anything or reaching for a separate device.

Why it matters:

Consistent settings and uninterrupted sessions lead to better outcomes. Control at your fingertips means you can fine-tune without breaking focus.

Session Timer Management

Set your session length (1–30 minutes) using the + / – buttons. The LCD display counts down in real time – auto-off when complete. No timers to set separately, no interrupted flow.

Result: Every session runs exactly as planned – start to finish, without manual intervention.

How to use

- Attach the wearable, then clip on the Battery Pack via the magnetic connectors
- Hold Power to turn on — LCD will illuminate
- From the home screen, use + / – to set your session timer (1–30 min)
- Press Mode and use + / – to select Red / NIR / Combined
- Press Mode again to set intensity, then again to set pulse
- Press Start — session begins immediately


RECOMMENDED USE

Session length	10–30 minutes per area
Frequency	Once daily for active recovery; or as directed by a health professional
Modes	Combined (Red + NIR) for joint recovery; Pulse for acute pain phases
Attachment	Secure magnetic fit — no readjustment needed during session

hello@flowlife.com
www.flowlife.com
Flowlife AB
Runebergsgatan 7
114 30 Stockholm
Sweden

Manual

Flowlight Wearable Battery pack

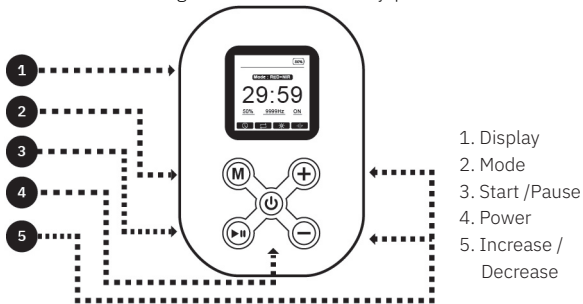Brochure

Flowlight Wearable Battery Pack

Flowlife

ENGLISH

INTRODUCTION

Welcome to your Flowlight Wearable Battery pack. This unit is the control centre and power source for your Flowlight Wearable recovery system. Read all instructions before first use and retain this manual for future reference.

SAFETY INSTRUCTIONS

Safety Hazards – Reduce the Risk of Electric Shock

- Use only the included USB-C cable or a certified USB-C power source for charging — do not use damaged cables or non-standard chargers
- Do not charge or operate the controller in wet conditions or near water
- Disconnect the USB-C cable immediately if the cable, port, or controller appears damaged
- Do not attempt to charge via any port other than the designated USB-C charging port
- Do not disassemble or attempt to service the internal charging circuitry
- Keep the USB-C port clean and dry — do not insert foreign objects
- Unplug from the power source when not charging

Danger – Reduce the Risk of Burns, Fire, or Injury

- Do not expose the Battery Pack to temperatures above 45°C or below -10°C during storage or use

- Do not place the controller on soft surfaces that may block ventilation during charging
- Do not crush, pierce, or place heavy pressure on the unit — risk of internal battery damage
- Stop use immediately if the unit becomes unusually hot, swells, or emits an unusual smell
- Do not use if the casing is cracked or the magnetic connectors appear damaged
- Keep away from flammable materials during charging
- Do not use near aerosol sprays, pure oxygen environments, or flammable gases
- Do not attempt self-repair — contact support@flowlife.com for servicing
- Keep out of reach of children under 15
- Do not leave charging unattended for extended periods in unsafe environments

Light Safety

- The Battery Pack itself does not emit therapeutic light — light is emitted by the wearable device it controls
- Follow all light safety instructions provided in the wearable device manual

General Safety

- Ensure the magnetic connectors are fully seated before starting a session
- Do not use the controller with non-Flowlife wearable devices
- Check that the LCD display activates correctly on power-on before each session
- If the display fails to activate or buttons are unresponsive, discontinue use and contact support

Risk Groups – Consult a Doctor Before Use

- Pregnant or nursing – follow guidance in the wearable device manual regarding light therapy use
- Active implanted electronic devices (pacemaker, ICD) – keep the magnetic controller away from implant sites; consult your doctor before use
- Any medical condition affecting the treatment area – consult your doctor before beginning therapy sessions
- Epilepsy or sensitivity to flashing light – exercise caution when using Pulse mode; consult a doctor before use
- Photosensitising medications – consult your doctor before use
- Under 15 years old – use only under supervision of a healthcare professional
- Any other health concern – consult your doctor before starting

INSTRUCTIONS FOR USE

Charging the Controller

1. Connect the USB-C cable to the charging port on the Battery Pack
2. Connect the other end to a USB-C wall adapter or power bank (not included)
3. LED indicator shows 100% when fully charged (approx. 2.5 hours from empty)
4. Disconnect after full charge

Attaching to the Wearable

1. Fit the wearable onto the target area and secure the straps
2. Align the magnetic connectors on the Battery Pack with the wearable connection points
3. The magnets will seat firmly – ensure a secure connection before powering on

Running a Session

1. Hold Power to turn on – LCD screen illuminates
2. Press Mode – use + / – to select mode: Red / NIR / Combined
3. Press Mode again – use + / – to set intensity level
4. Press Mode again – use + / – to set your Hz level
5. In the home screen view, press + / – to set session timer (1–30 min)
6. Press Start to begin – LCD counts down the session timer
7. Press Start/Pause at any time to pause or resume
8. Device powers off automatically when the session timer reaches zero

After Your Session

1. Hold Power to turn off if not already off
2. Detach the Battery Pack from the wearable by releasing the magnetic connectors
3. Charge the Battery Pack as needed before the next session

BATTERY & CHARGING

- Built-in rechargeable battery — 6000 mAh
- Charging input: USB-C (5V, 1.0A)
- Includes EU power adapter and USB-C cable
- Charging time: approx. 2.5 hours from empty to full
- LED/LCD indicator shows 100% when fully charged
- Do not use non-USB-C power sources for charging
- For storage over 1 week, keep battery at approx. 50–70% charge level
- Avoid full discharge before recharging when possible

CLEANING & STORAGE

- Wipe the controller exterior with a soft, dry cloth
- Do not expose to water, solvents, or cleaning sprays
- Keep the USB-C port and magnetic connectors free of debris — use a dry brush if needed
- Inspect the magnetic connectors before each session — clean contact surfaces gently if needed
- Store in a cool, dry location away from direct sunlight
- Do not store at temperatures below -10°C or above 45°C
- Do not stack heavy objects on top of the controller during storage

TECHNICAL SPECIFICATIONS

Specification	Value
Product Name	Flowlight Wearable Battery pack
Compatible with	Flowlight Wearable Knee 60 Two Waves, Flowlight Wearable Foot 60 Two Waves
Battery Capacity	6000 mAh
Charging Input	USB-C — 5V $\overline{=}$ 1.0A
Includes	EU power adapter + Type-C USB cable
Charging Time	Approx. 2.5 hours (empty to full)
Display	LCD screen — shows mode, intensity, and session timer
Controls	5 buttons: Power, Mode, Start/Pause, +, -
Session Timer	1–30 min (adjustable)
Light Modes	Red / NIR / Combined / Pulse (1–1000 Hz)
Pulse Frequency	1–1000 Hz (adjustable)
Connection	Magnetic clip — attaches directly to Flowlight wearable
Shell Material	ABS plastic

Dimensions	8.5 × 6 × 4 cm
Weight	0.2 kg
Certifications	CE, UKCA
Warranty	2 years

WARRANTY & SUPPORT

Flowlight Wearable Battery pack includes a 2-year warranty covering manufacturing defects from the date of purchase.

Warranty does not cover:

- Physical damage from misuse, drops, or accidents
- Water damage or use in wet environments
- Use with non-original or damaged USB-C cables
- Unauthorised modifications or repair attempts
- Battery degradation from improper charging habits
- Missing proof of purchase

For questions or service, contact: support@flowlife.com | www.flowlife.com

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