Flowlight Laser Mask Ultra

FLOWLIFE

Brochure

Flowlight Laser Mask Ultra

Red Light Therapy Mask

Table of Content

4 —	Manifest
6 —	Why Flowlight Laser Mask Ultra
7 —	Accessories
8 —	Flowlight Laser Mask Ultra - A Comprehensive Guide
10 —	Starting your Journey
16 —	Sleep Quality & Circadian Support
20 —	Stress Relief & Nervous System Reset
22 —	Energy & Cognitive Clarity
24 —	Skin Health & Facial Recovery
26 —	Contact Information
28 —	Manual
	• English • Svenska • Danska • Norsk • Suomi • Deutsch • Polski • Français • Italiano • Español • Nederlands

Manifest

In a world where limits are meant to be shattered, we fuel the fire of relentless pioneers. Our brand champions those who break the mold and embrace the extraordinary, symbolizing a life lived without compromise.

Committed to relentless innovation, we design for those who demand excellence and performance, understanding that the journey into the unknown often leads to unforgettable achievements.

We want to power dreams, aspirations, and triumphs of the audacious, sparking the passions of the unstoppable.


For those who redefine boundaries and explore new horizons, we are Flowlife. Made Possible.

Why Flowlight Laser Mask Ultra?

- **Precision Laser Therapy for Advanced Recovery** - Flowlight Laser Mask Ultra uses clinical-grade diodes at 660 nm, 850 nm, and 1060 nm—clinically validated wavelengths that target deeper tissue layers, enhance mitochondrial activity, and regulate inflammation. A 20° beam angle concentrates light into the face's vascular and nerve-dense zones, positioning the mask as an effective tool for full-body autonomic balance.
- **High-Density Laser Array with Deep Penetration** - With 164 laser-emitting diodes, the mask ensures full facial coverage. Each wavelength is placed for optimal surface and subdermal activation. Irradiance reaches ~90 mW/cm²—three times higher than standard LED masks—delivering more energy, deeper effects, and faster outcomes for muscle recovery, stress relief, and sleep support.
- **Structured Recovery, Customizable Intensity** - A wireless controller lets you choose between three light modes (Red + NIR, Red only, or NIR only), three timer options (5, 10, or 15 minutes), and three intensity levels. Whether used to support circadian rhythm, reduce neural tension, or recover from training, the protocol adapts to your needs.
- **Compact Power, Clinical Depth** - This is precision light therapy in wearable form. Flowlight Laser Mask Ultra complements full-panel systems by targeting cranial nerves, vagal points, and hormonal centers in the face. Designed for daily use, it delivers clinical-level outcomes in a compact format—without compromise.


Accessories

Flowlight Laser Mask Ultra comes fully equipped for high-performance light therapy wherever you are. Every component is selected for precision, portability, and safe operation—whether you're integrating it into a daily recovery ritual or using it as a high-efficacy complement to full-body panel therapy.

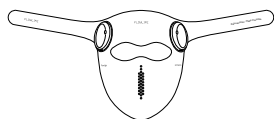Flowlight Laser Mask Ultra

The core device. Features focused-angle laser diodes (660 nm, 850 nm, 1060 nm) for deeper, more precise light penetration than LED. Can be activated via onboard button or controlled wirelessly. Ergonomic silicone design ensures full facial contact and targeted stimulation of cranial and parasympathetic zones.

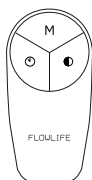Wireless Controller

Controls all sessions with three light modes, three durations (5, 10, 15 min), and three intensity levels. Tactile buttons and LED indicators make it easy to confirm settings. Real-time feedback ensures reliable pairing with the mask.


USB-C Charging Cable

Use to charge your Flowmask.

Protective Carry Pouch

Slim, durable, and structured for storage or travel. Keeps the mask and accessories protected without bulk.

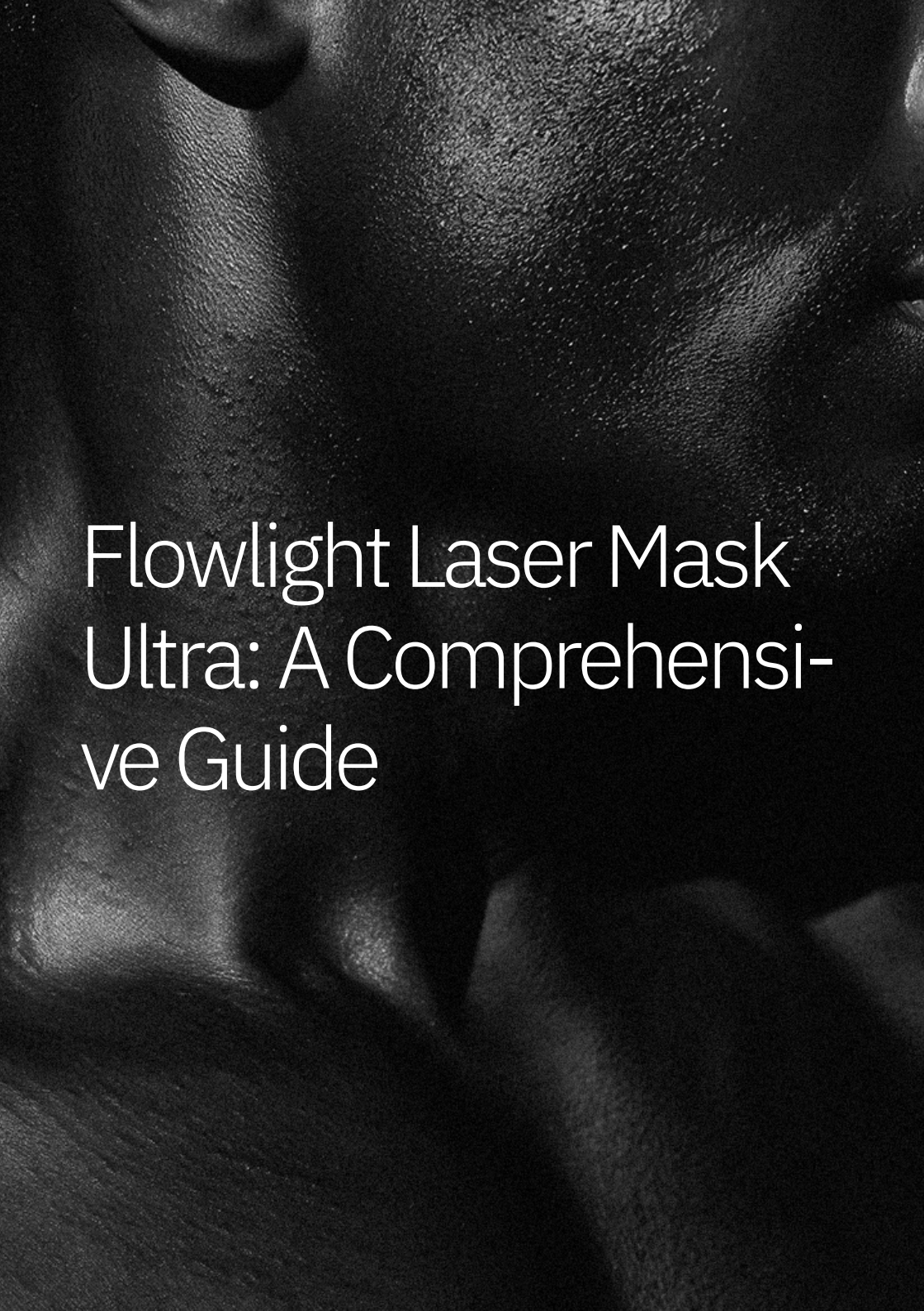A high-contrast, black and white close-up photograph of a person's face. The person is wearing a dark, textured laser mask that covers the entire face. The mask has a fine, pebbled texture. The lighting is dramatic, highlighting the contours of the face and the texture of the mask. The person's eyes are closed, and their mouth is slightly open, showing a small gap in the mask. The overall mood is mysterious and futuristic.

Flowlight Laser Mask Ultra: A Comprehensive Guide

To make the most of your sessions, it's essential to understand how Flowlight Laser Mask Ultra works and how to use it effectively, both with and without the wireless controller. This guide walks you through the entire experience, from activation to storage.

Preparing for Your Session

Start with clean, dry skin

Ensure your face is clean and free from make-up, oil, and lotion. Clean skin maximizes light absorption and minimizes interference with the mask's optical output.

No need for eye protection

Flowlight Laser Mask Ultra is designed with a safe, clinically aligned wavelength distribution. Most users prefer to close their eyes during treatment for comfort, but protective glasses are not required.

Starting Your Session

Activate Default Mode (manual use)

Hold the small button located on the top-right inner side (above your right eye when worn) for 3 seconds. This activates the device and the default mode.

- 15-minute treatment.
- Combined 660 nm red, 850 nm and 1060 nm near-infrared (NIR wavelengths).
- Intensity level 3 (highest).

Note: Without the controller, only the Default Mode is available. You cannot adjust mode, intensity, or time.

Understanding the Modes

To customize your session, connect the wireless controller. This unlocks complete control of the following:

- Mode: Red + NIR / Red only / NIR only.
- Timer: 5, 10, or 15 minutes.
- Intensity: Level 1 (Low), 2 (Medium), or 3 (High).

1. Red + NIR (Default Mode)

Wavelengths: 660 nm, 850 nm, 1060nm

Delivers broad-spectrum stimulation from superficial dermal layers to deeper neurological tissues. Supports systemic recovery, mitochondrial repair, and the regulation of inflammation.

2. Red Only (660 nm)

Optimized for skin-level applications. Stimulates collagen production, accelerates wound healing, and supports mood regulation.

3. Near-Infrared Only (850 nm, 1060nm)

Penetrates deeper into the facial nerves, vasculature, and connective tissues. Ideal for nervous system regulation, hormonal balance, and fatigue relief.

Intensity Levels

Each mode can be adjusted to 3 intensity levels via the controller. Higher levels increase irradiance but do not change wavelength.

Timer Settings

Available durations are 5, 10, or 15 minutes. The mask powers off automatically at the end of the session.

Visual & Sensory Mode Indicators

To help confirm which light mode is active, the Flowlight Laser Mask Ultra uses a built-in indicator light positioned on the outside of the mask's right-side module:

- White light = Mode 1 (Red + NIR)
- Solid red light = Mode 2 (Red only)
- Flashing red light = Mode 3 (NIR only)

In addition to the external indicator, the mask offers sensory and visual cues during treatment:

- Red-only mode emits a visible red light.
- Red + NIR mode emits the same red light, often perceived as slightly deeper due to the combined wavelengths.
- NIR-only mode may appear inactive to the eye, as near-infrared light is invisible. However, the mask will still emit gentle warmth and remain fully active. Use the indicator light to confirm mode selection.

Using the Wireless Controller

The Flowlight Laser Mask Ultra controller allows for complete customization of your treatment. It is wireless, intuitive, and syncs directly with the mask using both audible and visual feedback.

Initial Activation

After inserting the battery into the controller, hold the mode button "M" to activate the remote control.

Button Layout & Functions

Mode Button (top middle)

- Press once to cycle between light modes:
 - Mode 1: Red + NIR
 - Mode 2: Red only
 - Mode 3: NIR only
- The selected mode is confirmed by the indicator light on the outer surface of the mask's right-side module.

Note: This button does not power the mask on or off—it only changes the light mode.

Time Button (bottom left)

- Press to select treatment duration: 5, 10, or 15 minutes.
- The LED display shows 1, 2, or 3 lights, corresponding to the selected duration.

Intensity Button (bottom right)

- Press to set intensity level: Low (1), Medium (2), or High (3).
- The LED display shows 1, 2, or 3 lights to reflect the selected brightness.

Audible Feedback

Every adjustment—mode, timer, or intensity—is confirmed by a soft beep from the mask, ensuring real-time responsiveness even when the visual indicator is not in view.

During Your Session

Relax and disconnect

While Flowlight Laser Mask Ultra is built for portability, practical sessions require stillness. Recline, close your eyes, and allow the therapy to take effect without distraction. Consider pairing with breathwork or mindfulness for amplified effects.

Monitor comfort

Adjust settings via the controller if needed. Most users find Level 2 or 3 intensity ideal after initial acclimatization.

After Your Session

Automatic shutoff

The mask turns off automatically when the session ends. You can also stop it early by holding the power button on the mask.

Cleaning and storage

Gently wipe the device with a dry microfiber cloth. Store it in the provided carry pouch to keep it clean and safe.

Charging

Use the included USB-C charger to recharge the mask. Use two AAA-Battery for the remote control. One full charge typically supports multiple sessions.

Laser Mask Ultra

Sleep Quality & Circadian Support

Optimizing both ends of the sleep cycle—
morning activation and evening regulation.

Morning: Set Your Internal Clock

Why it matters

Your circadian rhythm is primarily regulated by light. Early exposure to red light tells your brain it's time to be awake, setting the internal clock for melatonin production 12–14 hours later. This improves your ability to fall asleep naturally in the evening.

How it works

- Red light (660 nm) suppresses melatonin and elevates natural cortisol to promote wakefulness.
- Enhances morning mood, energy, and mental clarity.
- Strengthens circadian timing for more stable nighttime regulation.

Recommended Settings

Mode

- Red Only or Combined Mode (Red + Near-Infrared).

Intensity

- Level 2–3.

Duration

- 10–15 minutes.

Time of Day

- Within 1 hour of waking.

Frequency

- Daily, especially in darker seasons or after travel.

Evening: Wind Down for Restorative Sleep

Why it matters

In the evening, your nervous system should transition from a state of alertness to one of relaxation. Near-infrared light supports this by calming facial tension and triggering parasympathetic activity—without stimulating the brain as bright blue light would.

How it works

- Near-infrared light (850 + 1060 nm) stimulates vagal and cranial nerve pathways.
- Reduces facial muscle tension and stress accumulation.
- Supports melatonin production by avoiding blue light suppression.

Recommended Settings

Mode

- Near-Infrared Only (for deep relaxation).

Intensity

- Start at Level 1–2; increase to Level 3 if well-tolerated.

Duration


- 10–15 minutes


Time of Day

- 30–60 minutes before bed.

Frequency

- Nightly or after high-stress evenings.

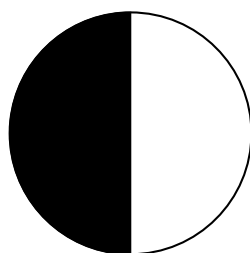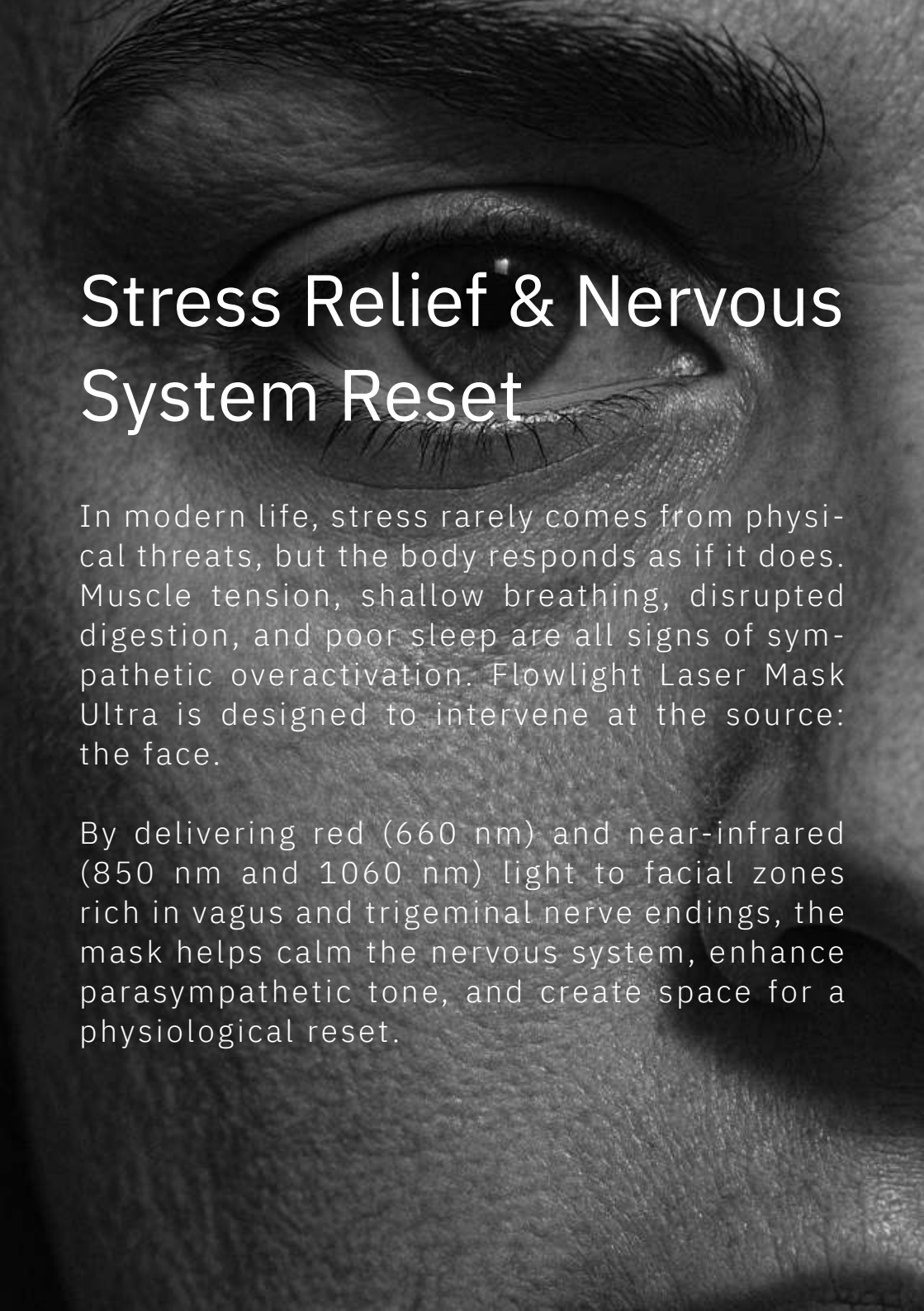Stress Relief & Nervous System Reset

In modern life, stress rarely comes from physical threats, but the body responds as if it does. Muscle tension, shallow breathing, disrupted digestion, and poor sleep are all signs of sympathetic overactivation. Flowlight Laser Mask Ultra is designed to intervene at the source: the face.

By delivering red (660 nm) and near-infrared (850 nm and 1060 nm) light to facial zones rich in vagus and trigeminal nerve endings, the mask helps calm the nervous system, enhance parasympathetic tone, and create space for a physiological reset.

How it works

- Light absorption through cranial nerve branches (especially supraorbital and infraorbital) activates parasympathetic pathways.
- This slows the heart rate, lowers blood pressure, and reduces the levels of circulating stress hormones.
- Regular exposure may improve vagal tone, heart rate variability, and emotional regulation.
- Many users report a sense of calm afterward, often accompanied by improved sleep or digestive ease.

Recommended Settings

Mode

- Red + NIR (660 + 850 + 1060 nm) for layered stimulation of both surface and deep nerve endings..

Intensity

- Level 1–2 (gentle input is more effective for nervous system regulation).

Duration

- 10–15 minutes

Time of Day

- Afternoon or early evening—avoid use during multitasking or high-activity.

Frequency

- 3–7 times per week, depending on baseline stress levels.

Energy & Cognitive Clarity

Brain fog, low motivation, and midday crashes often stem from disrupted circadian rhythms and sluggish mitochondrial function. The face—rich in cranial nerve endings and vascular density—serves as a gateway to systemic energy balance. Flowlight Laser Mask Ultra delivers precision-targeted red and near-infrared wavelengths to stimulate cellular energy, oxygenation, and cognitive performance.

How it works

- Red light (660 nm) activates mitochondrial function in superficial tissues and modulates wakefulness via hypothalamic signaling.
- Near-infrared light (850 nm and 1060 nm) penetrates more deeply into the vasculature and cranial structures, supporting stimulation of the frontal cortex and systemic ATP production.
- These effects help stabilize cortisol rhythms and elevate energy without nervous system overstimulation.
- Many users experience enhanced focus, clearer thinking, and steadier daytime energy.

Recommended Settings

Mode

- Red only or Red + NIR (660 + 850 + 1060 nm) for full-spectrum activation.

Intensity

- Level 2–3 (depending on sensitivity and desired stimulation)

Duration

- 10–15 minutes.

Time of Day

- Within 2 hours of waking or early afternoon, avoid use late in the day.

Frequency

- 3–6 times per week, ideally paired with a morning ritual or high-focus block.

Skin Health & Facial Recovery

Flowlight Laser Mask Ultra delivers medically validated wavelengths—660 nm red light and near-infrared light at 850 nm and 1060 nm—to support collagen synthesis, calm inflammation, and rebuild skin vitality from the inside out. Unlike topical solutions, these wavelengths penetrate through the epidermis and dermis, triggering deep biological repair.

How it works

- 660 nm red light stimulates fibroblasts to produce more collagen and elastin, which are crucial for improving texture, firmness, and skin tone.
- Near-infrared light at 850 nm and 1060 nm enhances blood flow and cellular metabolism, accelerating recovery and reducing oxidative stress in deeper layers.
- Anti-inflammatory effects calm reactive skin and promote healing from breakouts, redness, or fatigue.
- Repeated exposure strengthens the skin's barrier and improves long-term resilience.

Recommended Settings

Mode

- Red only (for surface-level repair) or Red + NIR (for full-spectrum dermal and subdermal support)

Intensity

- Level 2 for daily skin maintenance
- Level 3 for post-stress treatment or visible rejuvenation

Duration

- 10–15 minutes per session.

Time of Day

- In the morning or after cleansing, avoid use just before direct sun exposure.

Frequency


- 5–7 times per week for active recovery; 3–4 times per week for maintenance.

hello@flowlife.com
www.flowlife.com
Luntmakargatan 90
113 51 Stockholm

Manual

Flowlight Laser Mask Ultra

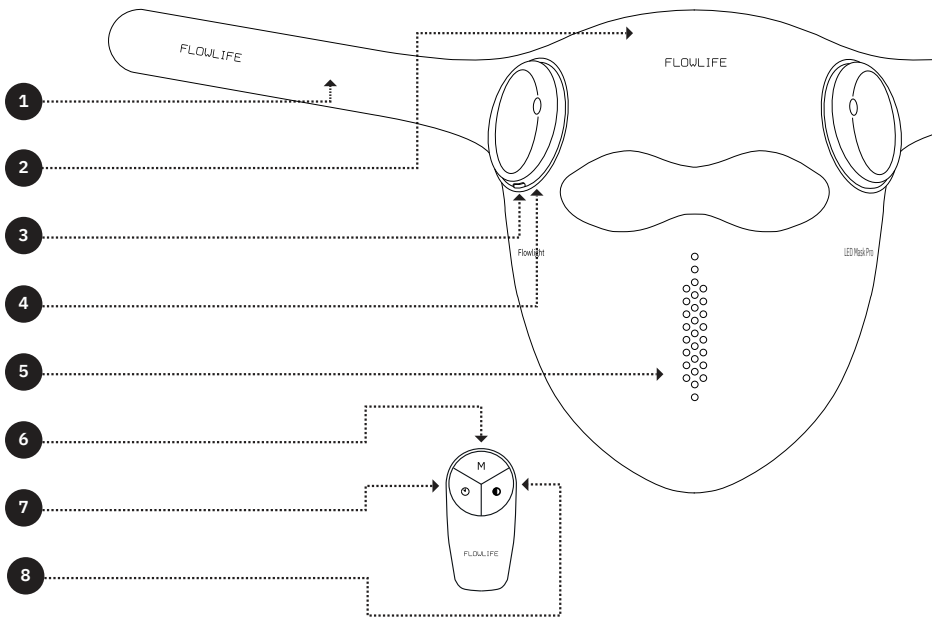1. Velcro Strap
2. Mask
3. On/Off
4. USB Charging Port
5. Breathing Vents
6. Mode Selection
7. Timer
8. Light Intensity

ENGLISH SAFETY INSTRUCTIONS

ELECTRICAL HAZARDS WARNINGS

- Never clean the mask while plugged in
- Do not use near water (bath, sink, etc.)
- Unplug immediately if it comes into contact with water

USAGE WARNINGS

- Never cover the mask during use
- Do not use on broken or irritated skin
- Only use original accessories
- Keep out of reach of children
- Never insert objects or attempt repair
- Store away from hot surfaces or damp areas
- Do not use while sleeping, driving, or operating machinery

STOP USE IF

- The product or cable is damaged
- It has been dropped or exposed to water

USAGE

- Ensure your face is clean and dry before use
- Place the mask securely on your face for full contact
- Fasten the Velcro strap behind your head to keep it in place during use
- Press and hold the small button on the inside top-right module (above your right eye when worn) for 3 seconds to activate the mask

DEFAULT MODE (MANUAL START)

- Duration: 15 minutes
- Wavelengths: 660 nm red + 850 nm & 1060 nm near-infrared
- Intensity: Level 3 (highest)

WIRELESS CONTROLLER (OPTIONAL FOR FULL CUSTOMIZATION)

- Insert 2 × AAA batteries into the wireless controller before first use
- After installing the batteries, all 6 indicator lights on the controller will flash 5 times—this confirms successful power-up and connectivity

MODE BUTTON (TOP MIDDLE)

- Press once to cycle between the three light modes:
- → Mode 1: Red + NIR
- → Mode 2: Red only
- → Mode 3: NIR only
- Mode indicator (on right-side module of mask):
- – White light = Mode 1 (Red + NIR)
- – Solid red = Mode 2 (Red only)
- – Flashing red = Mode 3 (NIR only)
- Audible beep confirms each selection
- Note: This button does not power the mask on or off

TIME BUTTON (BOTTOM LEFT)

- Press to set session duration: 5, 10, or 15 minutes
- LED indicators:
- – 1 light = 5 minutes
- – 2 lights = 10 minutes
- – 3 lights = 15 minutes

INTENSITY BUTTON (BOTTOM RIGHT)

- Press to set intensity level: Low, Medium, or High
- LED indicators:
- – 1 light = Low
- – 2 lights = Medium
- – 3 lights = High

DURING USE

- Sit or lie down comfortably
- Avoid multitasking
- Keep eyes closed if preferred

AFTER USE

- The mask will shut off automatically when the session ends
- To end manually, hold the button on the mask

CLEANING

- Wipe both inner and outer surfaces with a dry microfiber cloth after each use

CHARGING

- Use the included USB-C cable to charge the mask
- Charging indicator (on mask):
 - Red flashing light = Charging
 - White light = Fully charged

RISK GROUPS – ALWAYS CONSULT A DOCTOR IF:

- You have had recent burns, fever, infections, or photosensitivity
- You are pregnant (avoid the abdominal area)
- You are diagnosed with SLE, cancer, hyperthyroidism, or severe bleeding disorders
- You use photosensitizing medications
- You have eye disease or visual sensitivity
- You have a history of melasma or hyperpigmentation

DISCLAIMER

- This product is not a medical device
- It is designed for non-medical home and professional wellness applications
- Effects vary between individuals
- Consult a healthcare provider for specific concerns

ELECTRONIC WASTE & DISPOSAL

- This product contains a lithium-ion battery and electronic components
- Do not dispose of in household waste
- Recycle electronic waste in accordance with local regulations

WARRANTY

- A 2-year limited warranty, effective from the date of purchase, covers the Flowlight LED Mask Pro
- Covers defects in materials and workmanship under regular use
- Does not cover damage resulting from misuse, unauthorized repairs, accidents, or exposure to water
- Proof of purchase is required for warranty claims

TECHNICAL SPECIFICATIONS

- Wavelengths: 660nm (Red Laser), 850nm (NIR Laser), 1060nm (NIR Laser)
- Optical Power per Diode:
 - 660nm: 8mW
 - 850nm: 5mW
 - 1060nm: 10.5mW
- Beam Angle: 20°
- Irradiance on Skin: ~90 mW/cm²
- Battery: 1440mAh lithium-ion (40 min operation per charge)
- Number of Diodes: 164 (20 × 1060nm, 72 × 850nm, 72 × 660nm)

OPERATION MODES

- Mode 1: Red + NIR combined
- Mode 2: Red only
- Mode 3: NIR only

TIMER OPTIONS

- 5, 10, or 15 minutes

INTENSITY LEVELS

- Low, Medium, High

CHARGING INDICATOR LIGHT

- The indicator will flash red when charging. And full white light when fully charged.

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
USB-C OPLAADKABEL


- Nominale ingang: AC100–240V
- Nominale uitgang: DC 5V–1A

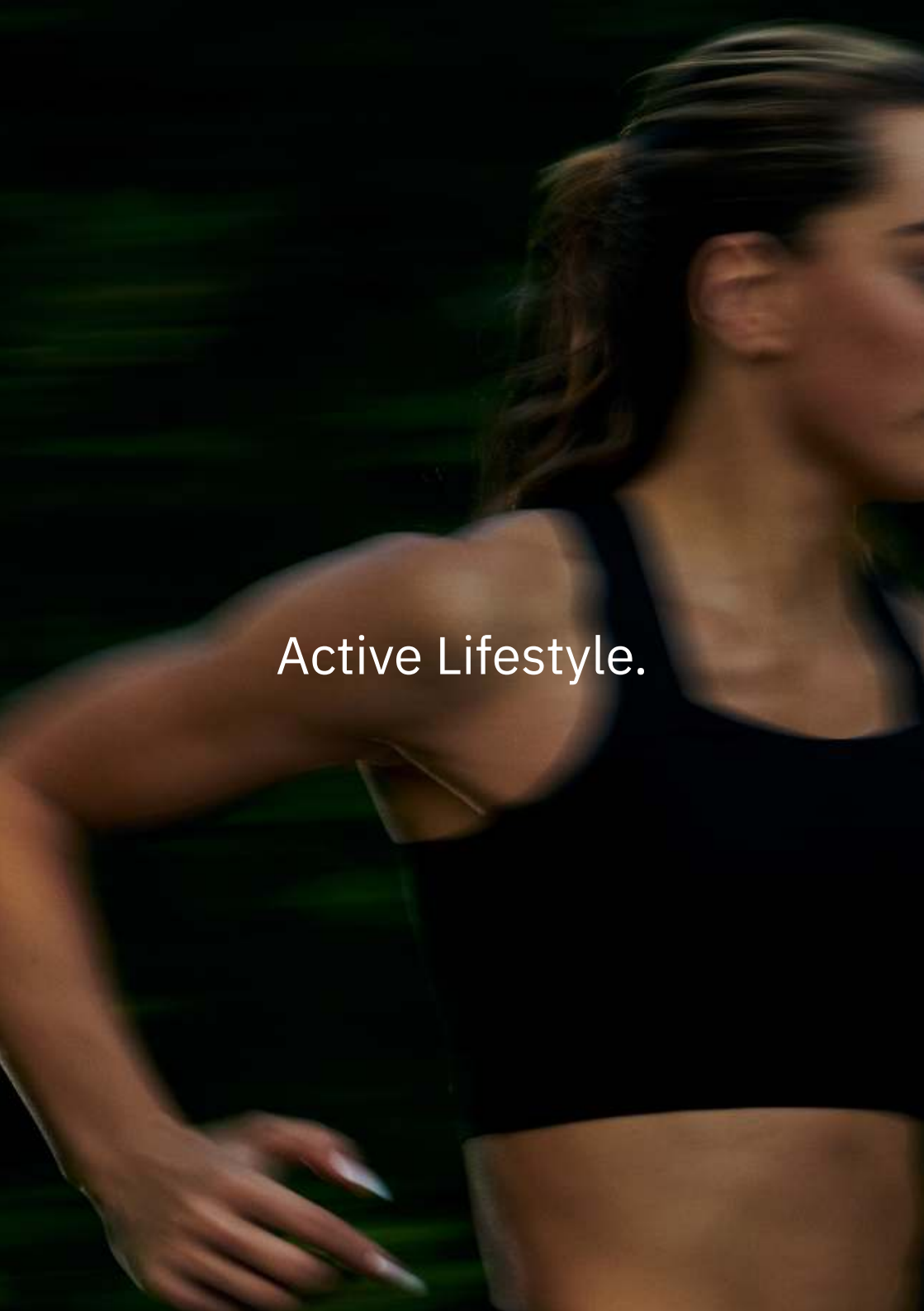A close-up, side-profile shot of a woman running at night. She is wearing a black athletic top. Her hair is dark and pulled back, with some strands blowing in the wind. The background is dark and blurred, suggesting a nighttime outdoor setting. The text "Active Lifestyle." is overlaid in white, sans-serif font in the center of the image.

Active Lifestyle.

Made Possible.

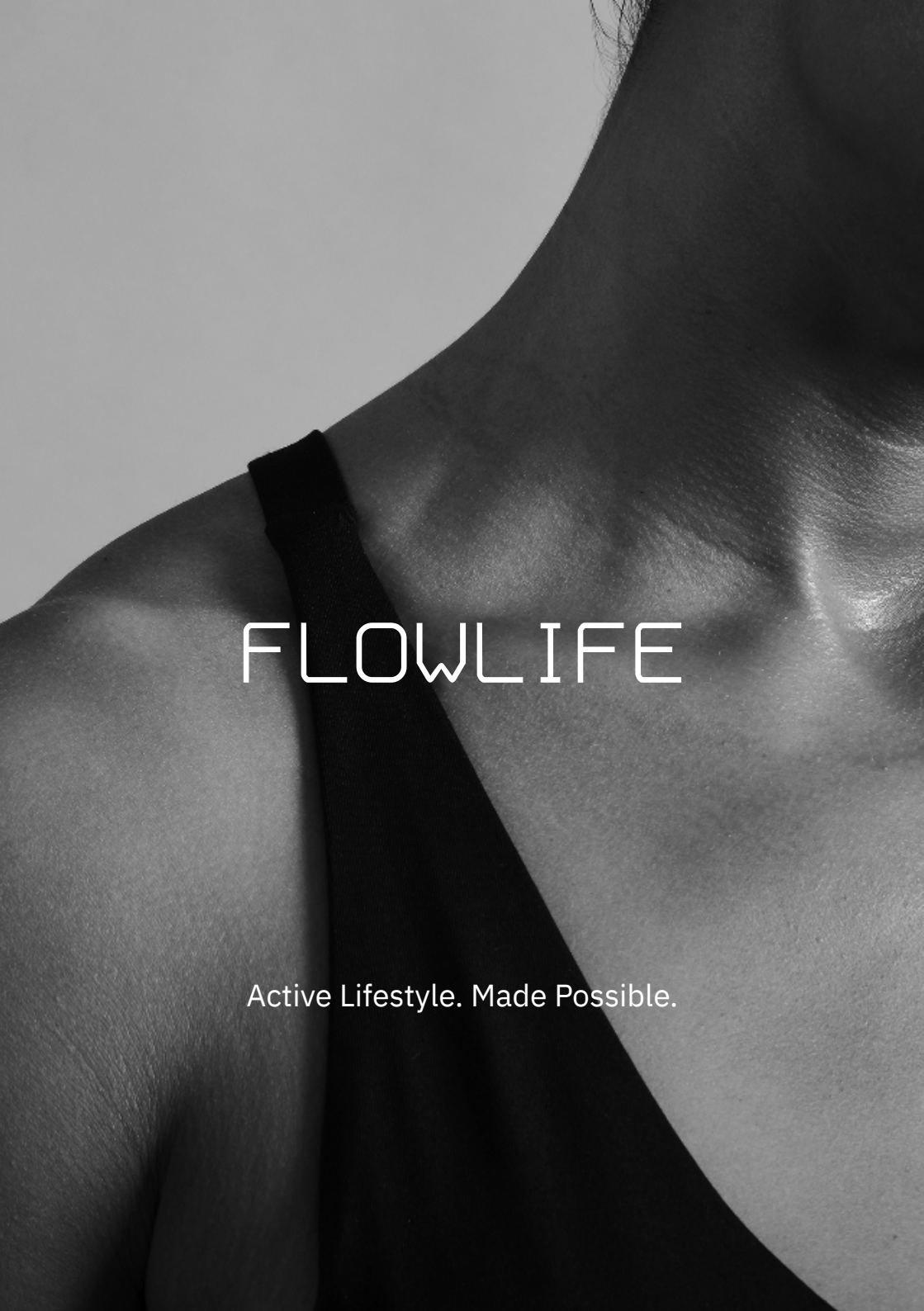A black and white photograph showing the back of a person's neck and shoulder. A black strap, likely from a sports bra or swimsuit, is visible across the shoulder. The skin is smooth and the lighting creates soft shadows, emphasizing the contours of the body.

FLOWLIFE

Active Lifestyle. Made Possible.