RAPTOR

SERVER ANALYSIS SOFTWARE

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TARTALOMJEGYZÉK

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Welcome to RAPTOR SERVER Analysis Software, the analysis software for the [ALEX multiplex test](#).
Please log in below to continue.

Log in with Azure AD / Office 365	Log in with Username & password
If your organization integrates the software into its Microsoft Azure/Office 365 user directory (single-sign-on) click the button below. You will be redirected to the Microsoft login screen.  Log in	Please enter your username and password. Your username Your password   Log in Forgot your password? No problem, you can request a new one.

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Log in with Username & password

Please enter your username and password.

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Contract approval required

Please review the following documents. You are required to confirm that you have read and agree with each document in order to use RAPTOR SERVER:

☐ **Bitte lesen** (last updated 23.11.2023 18:00)

☐ **RAPTOR SERVER IFU** (last updated 04.12.2023 11:05)

☐ **Terms of Service** (last updated 22.11.2023 01:00)

Confirm

By clicking the button, I confirm that I have reviewed these documents and agree.

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
New Measurements					Last 7 days	
QC	Product type / Source	QR-Code	Sample Code	Analysis Date		
Failed	Research MAX 45k	99AFG0F7	NK FOX	30/03/2023 20:50	Details	Report
Failed	Research MAX 45k	99AFG0F8	N.G.	30/03/2023 20:50	Details	Report
Failed	Research MAX 45k	99AFG0F9	E18735	30/03/2023 20:50	Details	Report
Failed	Research MAX 45k	99AFG0FA	20220325w56	30/03/2023 20:50	Details	Report
Failed	Research MAX 45k	99AFG0EE	HM100038	30/03/2023 20:50	Details	Report
Failed	Research MAX 45k	99AFG04D	786118	30/03/2023 20:37	Details	Report
Failed	Research MAX 45k	99AFG016	HM100225	30/03/2023 20:37	Details	Report
Failed	Research MAX 45k	99AFG0CA	E17685_M	30/03/2023 20:37	Details	Report
Failed	Research MAX 45k	99AFG0CC	E17685_M	30/03/2023 20:37	Details	Report
Failed	Research MAX 45k	99AFG015	E17685_F	30/03/2023 20:37	Details	Report
In total there are 270 new measurements in the selected time period.						
View all new measurements						

Approved Measurements						
QR-Code	Product type / Source	Sample Code	Analysis Date			
02AVJ227	ALEX2	Anna LMacey	28/03/2023 21:17	Details	Report	
02ARA379	ALEX2 IX	SM_28032023	28/03/2023 14:29	Details	Report	
02ARA393	ALEX2 IX	RG_28032023	28/03/2023 14:26	Details	Report	
02AVA31F	ALEX2 MAX 45k	PR_Mar2023_11	21/03/2023 15:05	Details	Report	
02AVI3B9	ALEX2 MAX 45k	KF	17/03/2023 14:46	Details	Report	
02AYT337	ALEX2 MAX 45k	2B46B46+CCD	14/03/2023 19:06	Details	Report	
02ATL349	ALEX2 IX		22/02/2023 14:21	Details	Report	
02AVE383	ALEX2 MAX 45k	AO	03/02/2023 09:41	Details	Report	
02AXI29A	ALEX2 MAX 1k	OW	19/01/2023 19:12	Details	Report	
02AUY3AE	ALEX2 MAX 45k	BSZ_Okt2022_8	18/10/2022 18:12	Details	Report	
View all approved measurements						

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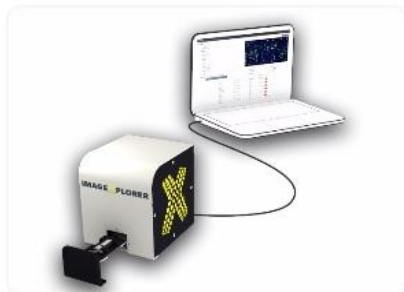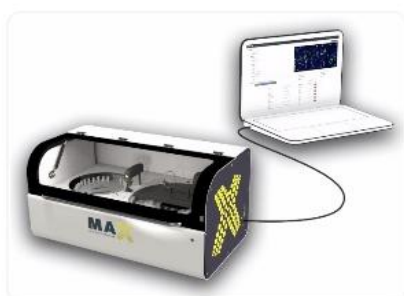
9K-Test-Instrument

Not initialized

[Initialize](#) [Shutdown](#) [Monthly Maintenance](#)

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E-Mail address:

Language: English (en)

The language is used to determine the display language.

Culture: English (United States) (en)

The culture is used to format numbers and dates.

Timezone: (UTC+01:00) Amsterdam, Berlin, Bern, Rome, Stockholm, Vienna

Configure your timezone to see date and times in your local time instead of UTC.

Multi-Factor Authentication

[Manage MFA](#)

Use multi-factor authentication (MFA) to significantly increase security. It is highly recommended to enable MFA; this adds a second layer of protection against phishing attacks or weak passwords.


Configuration

Research spots: Show research spots

If enabled, "research use only" arrays. This does not affect the

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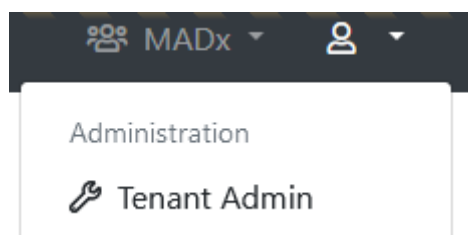General configuration

*** Name:** MADx

Default timezone: (UTC+01:00) Amsterdam, Berlin, Bern, R...
Configure a timezone to see date and times in (instead of UTC). Please note that users can choose to use a different timezone.

Default report template (ALEX2): ALEX2

Default report template (ALEX3, ALEX Air, ALEX Food): check 1.21

Default report template (FOX): FOX (2025)

Default report template (PAX): PAX

Default report language: English (en)

Default ALEX MOD: < No MOD >

Default FOX MOD: < No MOD >

Default PAX MOD: < No MOD >

Force MOD: False
Requires a measurement to always have an MOD assigned. If this setting is enabled, a measurement cannot be approved without an MOD assigned.

Hide Global MOD: False
If this setting is enabled, only the tenant-specific MOD will be available to select, and the global system-wide MOD will be hidden.

Display MOD Excluded Spots As Research Use Only: False
If this setting is enabled, spots excluded by an MOD will not be hidden from the measurement image, but displayed as a "RUO spot".

MFA: ☐ Force MFA
If this setting is enabled, users are forced to set up an MFA on sign-in, in case they haven't yet.

MAX device default settings

*** Default Minimum RunId:** 10
When a new MAX device is created, this value will be used for its **Minimum Run Id**

Tenant-specific batch releases

Restrict this tenant to only allow the following batches to be scanned or to restrict certain batches. Enter one per line.

Use:

☒ Whitelist
☐ Blacklist

Whitelisted Batches:

Enter one batch code including product identifier per line, for example "02AAA"

Blacklisted Batches:

Enter one batch code including product identifier per line, for example "02AAA"

Company Details

*** Legal Entity Name:** MacroArray Diagnostics

*** Street and Number:** Street 4

*** Zip Code:** 1230

*** City:** New York

*** Country:** United States of America

*** Phone:** +4318652573

*** E-Mail:** office@macroarraydx.com

*** Legal Representative:** Dr. Christian Harwanegg

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You requested a new measurement from ImageXplorer "IX Trainingsraum". Please wait until the measurement is completed ...

Invalid batch detected! The measurement can't be processed because its batch is blocked.

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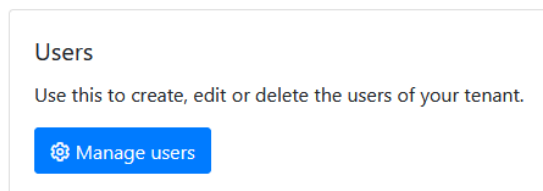
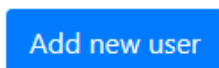
Add New User

General Information

Name:	New user	E-Mail address:	name@server.com
Status:	Enabled (user can log in)		

Login Credentials

Type:	Username and Password	Please select if the user should be validated against Azure AD (Office 365), or if you want to specify username and password yourself. It is recommended to use Azure AD whenever possible. This allows for increased security by supporting two-factor authentication. Also, it is more convenient for the user, because he or she does not have to remember a new username and password (single-sign-on).
Username:		
It is a good idea to use the users e-mail-address as his username, because the username needs to be unique across the entire system.		Password: Automatically send invitation e-mail

Access Rights

Preconfigured groups

Use these pre-configured groups to automatically assign a set of access rights to this user. You can still add or remove specific access rights.

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Tenants

This user is connected to the following tenants and has the following access rights:

☐ Tenants	Can share measurements	Can edit measurements	Can edit personal data	Can approve measurements	Can manage basic tenant configuration	Can manage API keys	Can manage ImageExplorers	Can view agent connections	Can manage users	Can manage users of Result Sharing Portal	Can manage report templates	Can manage measurement data completion	Can unapprove tenant specific batch releases	Can use a custom calibration string	Can execute MAX d	Can skip t
☐ Another MADx Tenant	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
☐ MADx Tenant	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓

Additionally this user can be added to the following tenants by checking the tenant name and hitting save afterwards. User rights w applied from current open tenant user.

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Additionally this user can be added to the following tenants by checking the tenant name and hitting save afterwards. User right applied from current open tenant user.

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Tenants

This user is connected to the following tenants and has the following access rights:

⚙️ Edit selected tenants rights

☐	Tenants	Can share measurements	Can edit measurements	Can edit personal data	Can approve measurements	Can manage basic tenant configuration	Can manage API keys	Can manage ImageXplorers
☒	Another MADx Tenant	✓	✓	✓	✓	✓	✓	✓

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Add new ImageXplorer

General information

Name:	
Key:	Cxltclg8okSwQ6CqJ4EZLWGnw2PdJ3
	This key is required when installing an Agent.
ImageXplorer Serialnumber:	
	Serialnumber of connected ImageXplorer.
Notes:	
Internal Notes:	
	For internal use only.
Disabled:	Not disabled

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Analysis Log

```
'PreProcessedPixels' took 128ms.
Using definition '02_ALEXII_IVD_v1.6.0'.
'ApplyArrayDefinition' took 74ms.
Parsing calibration string ...
'ApplyDefinitions' took 75ms.
'FindGrid' took 1ms.
'AlignFeatures' took 862ms.
'CropImage' took 1ms.
'QuantifyFeatures' took 231ms.
'SubtractBackground' took 681ms.
Parsing calibration string ...
Using calibrators from QR-Code.
Total IgG (tIgG) normalized by calculating 8,330999755859375E-07x*+0,05697999954223633x+-30,69269921875
(x=8515), which results in 515.
'ReQuantify' took 57ms.
```

QC Log

These are the values used for QC (ALEX v2): Curve Correction=0,9997370867894242, tIgG=8515, GD1=28085 (SNR=105.38), GD2=29369 (SNR=105.28), GD3=28569 (SNR=147.53), GD4=28305 (SNR=128.1). Reference values are as follows: Guide Dots must have a value above 5000; each Guide Dot spot must exist and be unique; Guide Dot spot must not be set to bad. Average background must be between 4000 and 16000. QC (ALEX v2) succeeded.

Background Log

Calculated background statistics for 164 pixels (from a total of 182, skipped smallest 9, and highest 9 pixels). The average is 9.419, the minimum 7973 and the maximum 11792 (minimum and maximum is calculated including skipped pixels). This results in an average difference of 628, which is an SD of 49,17 and a CV of 0,52%. For comparison, there is an overall total of 282944 pixels. The overall average is 14.071, the minimum 6.012 and the maximum 64.965.

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
Tissue & Source		Feature	kUA/L
Pollen Inhalation		E Aca m	≤ 0.10
Mites		E Aca s	7.11

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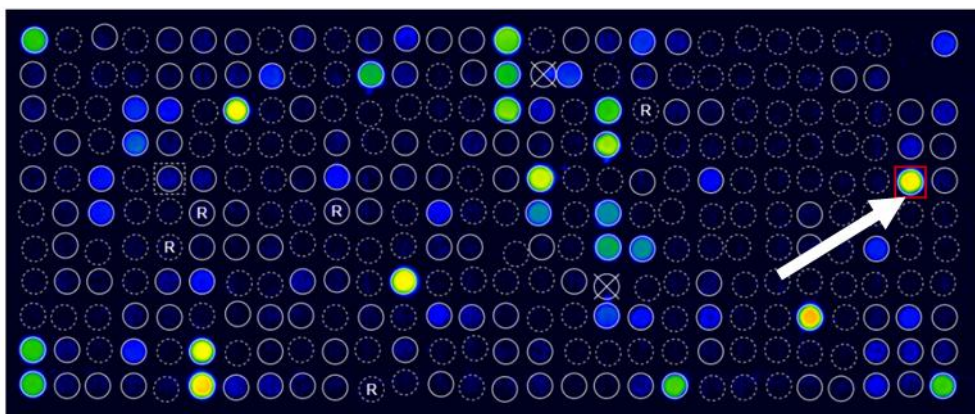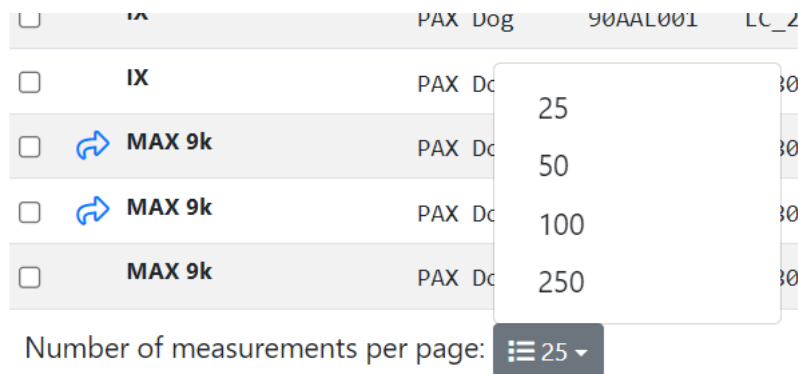
New Measurements

Search by QR-Code, patient, sample, tag, analysis date, approval date ...

Selected measurements

☐	Source	QC	Tags	Product type
☐	IX	✓ OK		ALEX2
☒	IX	✓ OK		ALEX2
☒	IX	✓ OK		ALEX2
☒	IX	✓ OK		ALEX2

Download CSV format
Download Excel format
Approve
Update interpretation
Choose Reporttemplate
Modify Tags
Move
Delete

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Re-process measurements

You are about to re-process the following measurements. undone, the existing analysis result will be irreversibly over

QR-Code	Sample Code
02APK3C8	Dummy8
02APK258	Dummy1
02APK07F	Dummy13

Re-process these measurements
Cancel

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New Measurements

Search by QR-Code, patient, sample, tag, analysis date, approval date ...

	Source	QC	Tags	Product type
☐	Recaptured	✓ OK	recaptured	ALEX1
☐	⚠ Recaptured	✓ OK	recaptured	ALEX1
☐	QR-Code was scanned more than once.		recaptured	ALEX1

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RAPTOR New measurements Approved measur

New Measurements

Search by QR-Code, patient, sample, tag, analysis date, appr

Selected measurements ▾

	Source	QC	Tags	Product type	QR-Code
☒	IX	✓ OK		ALEX2	02AMC180
☒	IX	✓ OK		ALEX2	02AMC078
☒	IX	✓ OK		ALEX2	02AKV2D9

Selected measurements ▾

- Download CSV format
- Download Excel format
- Approve**
- Update interpretation
- Re-Process
- » Move
- Delete

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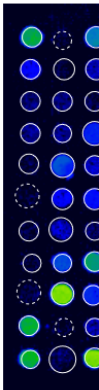
Approve measurements

You are about to approve the following measurements.

There are measurements with one or more warnings. Please make sure to double-check these measurements before approval.

QR-Code	Sample Code	Analysis Date
02AKV035		28.04.2022 11:14

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RAPTOR
New measurements
Approved measurements
Start measurement

Approved Measurements

Search by QR-Code, patient, sample, tag, analysis date ...

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RAPTOR
New measurements
Approved measurements
Start measurement
Automats
MADx Testtenant

Exports

Certain types of data exports (for example, the export of many PDF reports at once) take a few minutes to prepare for download. Here download your requested exports. Please note that you can keep working, all exports are done in the background and you can download them when they are completed.

Export requested	# of measurements included	Export completed
8/16/2021 11:07 AM	4	0 / 4
8/16/2021 11:06 AM	4	8/16/2021 11:06 AM
8/16/2021 11:04 AM	3	8/16/2021 11:04 AM

Download
Download
Download

Administration
Tenant Admin
Exports
Personal Settings
Log out
Support
Imprint
Download
Download

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Export requested	# of measurements included
06.07.2022 08:52	3

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Export multiple measurements

Please note that **Report PDF** and **ASTM** export is only a

QR-Code/Product/Batch		Sa
02ABC		
Analysis Date (from)	Analysis Date (to)	
TT . MM . JJJJ	TT . MM . JJJJ	
Tag	Status	
	All	
☐ Export CSV	☐ Export Excel	
Search Measurements		Cancel

Please note that **Report PDF** and **ASTM** export is only

QR-Code/Product/Batch

02aaa

02aaa

TT.MM.JJJJ

TT.MM.JJJ.

All

Export Excel

← Back

66 Measurements found that match your criteria

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Edit report template

Your changes have been saved successfully.

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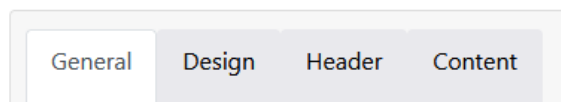
Configure report template "New report template"

General Design Header Content

Name: New report template This is the name of the template which shows up in the dropdown list while generating a report.

Base Template: Version 2 (ALEX2 only) A base template which specifies basic functionality must be chosen.

Preferred language: English (United States) (en-us) A default language can be specified which will be used if there isn't yet a language stored to a measurement's report.

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PDF naming

Setting Type: Manual

The filename of the report PDF can be customized. Either built automatically based on your selection below or manually adjusted in the order of your needs.

Name Pattern: Report_{PatientId}_{MeasurementQrCode}

Available Elements: {PatientId} ; {PatientName} ; {SampleId} ; {MeasurementQrCode} ; {DateAndTime}

--

PDF naming

Setting Type: Automatic

The filename of the report PDF can be customized. Either built automatically based on your selection below or manually adjusted in the order of your needs.

Include Patient ID: ☒

Include Patient Name: ☒

Include Sample ID: ☐

Include Measurement QR-Code: ☐

Include Date & Time: ☒

Footer

Footer text:

Will be printed in the footer if provided and rendered as HTML.

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Configure report template "New report template"

The screenshot shows the 'Design' tab of the 'Configure report template' interface. It includes the following elements:

- Tabs:** General, Design (active), Header, Content.
- PDF format:** Radio buttons for A4 (selected) and Letter. A description states: 'Specifies the paper size of the output PDF.'
- Font:** A dropdown menu showing 'Arial'. A description states: 'The font family of the whole report texts can be adjusted.'
- Background image:** A text input field and a 'Browse' button. A description states: 'If set, this image is used as background on the first page of the PDF.'
- Footer:** 'Maximum file size 1 MB'

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Third level headline color:

Logos

Logo (left): Browse
Maximum file size 1 MB
This logo will appear on the **top left corner** of the first page of the report.

Text (center):
This text will appear either on the left, right, or center of the top of the first page (depending on your configured logos). Supports the placeholders `{{ PATIENT_PORTAL_QRCODE }}` and `{{ PATIENT_PORTAL_URL }}` to dynamically insert measurement information.

Logo (right): Browse
Maximum file size 1 MB
This logo will appear on the **top right corner** of the first page of the report.

Small logo: Browse
Maximum file size 1 MB
This logo is used on all detail pages. If not set, the default logo is used.

Signatures

Left: Browse
Maximum file size 1 MB
This signature or image will appear on the last page, left aligned.

Center: Browse
Maximum file size 1 MB
This signature or image will appear on the last page, center aligned.

Right: Browse
Maximum file size 1 MB
This signature or image will appear on the last page, right aligned.

Indicators

Custom intensity indicator:
A custom SVG for displaying allergen intensity can be pasted here. Will be styled according to specified colors.

⚠ The SVG element must contain the attributes **width** and **height** with a minimum dimension of 16 for correct display in PDF.

🌿 SVGs paths which should be colored must require the **colorize** class as well as a fill defined in the **defs** segment of the SVG.

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General Design Header **Content**

Header content:

Left half		Right half	
Custom headline		Custom headline	
Patient Name	☐	Doctor Information	☐
Patient Name	☒	None	☐
Custom headline			
Date of Birth	☐	None	☐
Sample ID	☐	None	☐
QR-Code	☐	Additional Info text	☐
Analysis Date	☐	None	☐
Amount of tested allergens	☐	None	☐
Test method	☐	Quality Control text	☐

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General	Design	Header	Content
Allergens to include:		All allergens	Specifies which kind of allergens should be included in the report.
Allergens sort order:		Alphabetically by allergen	Specifies the sort order of the allergens.
Include Overview:		☒	Whether the allergen overview should be printed or not.
Include product type in overview:		☒	Whether the product type should be printed in the header of allergen overview.
Include Total IgE:		On the front page	Whether the information of total IgE should be included in the report and where.
Cross-reactive information:		Include text	This section will be displayed underneath the allergen list itself.
Include Interpretation:		Include manually	It can be specified if the Raven report interpretation information should be included in the report or not. This section is printed underneath the allergens overview.
i Additional information (below the list of results) Those informations are integrated into the report in an extra section underneath the allergen list.			
Approval date:		☐	
Assay date:		☐	
Doctor info:		☐	The information about the doctor which accompanied the measurement/patient.
Print date:		☒	When the report was generated.
Report info:		☐	The „additional info“ field belonging to the measurement.
Sample date:		☒	

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Configure report template "New report template - Duplicate"

General	Design	Header	Content
Name: New report template - Duplicate			
Support products: ALEX ²			
Preferred language: Deutsch (de)			

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*** Name:** MADx

Default timezone: (UTC+01:00) Amsterdam, Berlin, Bern, Rome, Stockholm, Vienna

Configure a timezone to see date and times in (instead of UTC). Please note that users can choose to use a different timezone.

Default report template (ALEX 1): Version 1 (ALEX1 & ALEX2)

Default report template (ALEX 2): Version 2 (ALEX2 only)

Default report language: English (United Kingdom) (en-gb)

Default MOD: < No MOD >

Require ImageXplorer testarray scan: Warn after 60 days, enforce after 90 days.

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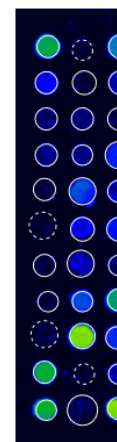
Template: Version 2 (ALEX2 only)

Language: English (en)

Usually this tenant's default template and language will be used when sharing or exporting measurements. You may change the template and language associated with this measurement by saving your current configuration.

You can download this report in PDF format once the measurement has been approved.

You have access to manage the report templates and settings via the tenant administration. Please note that any changes will affect all users.

MADX
MACROARRAY DIAGNOSTICS

MADX
MACROARRAY DIAGNOSTICS

PATIENT ID
PATIENT NAME
DATE OF BIRTH
SAMPLE ID
ASSAY CODE
02APH008
ANALYZED ON
11/04/2022
TESTED ALLERGENS
202
TEST METHOD CODE
ALEX2

REFERRING PHYSICIAN

NOTE
The internal QC (Fluorescence check for QC) was within acceptance range.

Lab report: Summary on detectable sensitisations

POLLEN
Grass Pollen

MICROORGANISMS
Fungal Spores & Yeast

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The screenshot shows the 'RAPTOR' web application interface. At the top, there are navigation links: 'New measurements', 'Approved measurements', and 'Start m...'. Below this is a section titled 'Approved Measurements' with a search bar labeled 'Search by QR-Code, patient, sample, tag, analysis date ...'. A dropdown menu 'Download selected measurements' is open, showing options: 'Report (PDF format)' (highlighted with a yellow circle), 'CSV format', 'Excel format', and 'ASTM format'. Below the menu is a table with columns: Source, Tags, Product type, QR-Code, and Sample C. The table contains three rows of data, all with 'IX' as the source and 'ALEX2' as the product type. The QR-Code for all rows is '02AKV39D'.

	Source	Tags	Product type	QR-Code	Sample C
☒	IX		ALEX2	02AKV39D	
☒	IX		ALEX2	02AKV39D	
☒	IX		ALEX2	02AKV39D	

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Tags:

MOD:

Report Template: Version 2 (ALEX2 only) - Deutsch (de)

[New measurements](#)
[Approved measurements](#)
[Start measurement](#)
[Automats](#)

Report

Template: Default

Language: English

[Save](#)
[Use defaults](#)

Usually this tenants default template and language will be used when sharing or exporting measurements. You may change the template and language associated with this measurement by saving your current configuration.

[Download PDF](#)

[Manage report template and settings](#)

You have access to manage the report templates and settings via the tenant administration. Please note that any changes will affect all users.

Patient ID:	
Patient Name:	
Date of Birth:	
Sample ID:	
Barcode:	02AAE
Sampled on:	
Tested on:	12/18/
Printed on:	12/19/
Note:	The internal QC (Plausibility)

Lab report: Summary on de


Pollen	Grass Pollen
	Tree Pollen
	Weed Pollen
Mites	House Dust M
Microorganisms	Fungal Spores
	Legumes
Plant-Based Food	Grain
	Spices
	Fruits
	Vegetables
	Nuts & Seeds
	Milk

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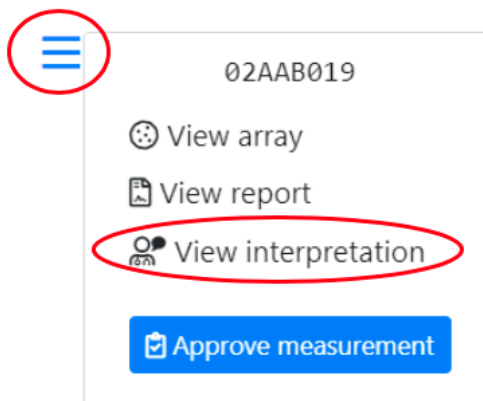Rule set

RAVEN Ruleset A for ALEX

A rule set defines the RAVEN or RAVEN² engine that is used to interpret results.

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AIT

Medication

MADx data declaration for smartRaven

☒ I accept to enter clinical data for using the full interpretation function of smartRAVEN

Patient has reported any symptoms

☐ Yes ☐ No

Additional information about the patient

☐ Atopic Dermatitis

☐ Allergen-specific immunotherapy (AIT) - ongoing, completed or discontinued

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[Request new interpretation](#)

This measurement already has an interpretation stored. Usually, it should not be necessary to request those again, but it could be necessary when the patients information has changed.

Interpretation

Raven Interpretation Summary

Sample Information

The sample was tested on ALEX² Barcode 02AAB019, interpretation date 12/19/2019.

Of the tested 295 allergens, 25 were above the detection limit of 0.3 kU_L/L. A sensitization can be an indicator of an allergy. For selected allergens, comments for interpretation are listed below.

Total IgE: 512 kU/L

The measured total IgE was 512 kU/L. Above a total IgE of 100 kU/L, allergy is likely.

Cross-Reactive Allergen Sensitization detected

Sensitization against allergen molecules which are markers of (broad) cross-reactivity between different allergen sources.

Detected cross-reactive allergen sensitizations:

- PR-10s: Aln g 1, Bet v 1, Cor a 1.0103, Fag s 1
- Polcalcins: Aln g 4, Phi p 7

PR-10 Proteins (PR10)

PR-10 inhalative: The major birch pollen allergen, Bet v 1, represents the prototype of all PR-10 allergens and is the primary sensitizer in regions with birch pollen exposure. The presence of PR-10 allergens in Fagales tree pollen explains IgE cross-reactivity between pollen from hazel, alder, beech, oak and hornbeam. PR-10 nutritive: PR-10 allergens in raw fruits, nuts, vegetable and legumes can induce oral allergy syndrome and sometimes severe allergic reactions in sensitized individuals. PR-10 allergens are not stable to heat and digestion.

Polcalcins (Ca++ binding proteins, PC)

Members of the PC family can cause inhalative symptoms. PC allergens can be found in pollen from grasses, trees and weeds. Inhalative symptoms manifest as allergic rhino-conjunctivitis and/or allergic asthma.

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The screenshot displays the 'Details' form in the MacroArray software. The form includes fields for Barcode (024E1852), Sample Code, Sampling Date (tt.mm.jjj), Assay Date (07.10.2020), Patient Code, Patient Name, Date of Birth (tt.mm.jjj), Info, Doctor Info, Tags, and MODs (Pollen allergens). The MODs field is circled in red. Below the form are buttons for Report, Save, and Save and approve. The bottom section shows the resulting spot array image, which is a grid of spots with blue and green colors. The image is labeled 'Previous' and 'Next' at the top right. Below the image, the Received date and time (07.10.2020 14:44) and the Analysis date and time (07.10.2020 14:44) are displayed. The ImageXplorer ID (23156891) and the Automatic QC status (OK) are also shown.

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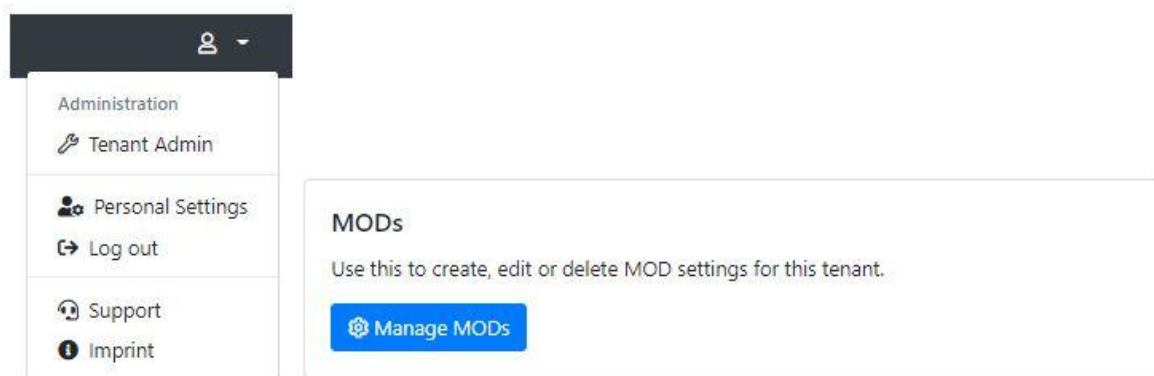
Type:

Deactivated:

You cannot delete an MOD once it is in use. But you can deactivate it, so that it's not available for new measurements anymore.

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*** Name:**

Default timezone:

Configure a timezone to see date and times in (instead of UTC). Please note that users can choose to use a different timezone.

Default report template (ALEX 1):

Default report template (ALEX 2):

Default report template (FOX):

Default report template (PAX):

Default report language:

Default ALEX2 MOD:

Default FOX MOD:

Default PAX MOD:

Force MOD:

Requires a measurement to always have an MOD assigned. If this setting is enabled, a measurement cannot be approved without an MOD assigned.

Hide Global MOD:

If this setting is enabled, only the tenant-specific MOD will be available to select, and the global system-wide MOD will be hidden.

Display MOD Excluded Spots As Research Use Only:

If this setting is enabled, spots excluded by an MOD will not be hidden from the measurement image, but displayed as a "RUO spot".

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MOD:

Report Template:

Force MOD:
Requires a measurement to always have an MOD assigned. If this setting is enabled, a measurement cannot be approved without an MOD assigned.

Hide Global MOD:
If this setting is enabled, only the tenant-specific MOD will be available to select, and the global system-wide MOD will be hidden.

Display MOD Excluded Spots As Research Use Only:
If this setting is enabled, spots excluded by an MOD will not be hidden from the measurement image, but displayed as a "RUO spot".

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Force MOD:

False

Requires a measurement to always have an MOD assigned. If this setting is enabled, a measurement cannot be approved without an MOD assigned.

Hide Global MOD:

False

If this setting is enabled, only the tenant-specific MOD will be available to select, and the global system-wide MOD will be hidden.

Display MOD Excluded Spots As Research Use Only:

True

If this setting is enabled, spots excluded by an MOD will not be hidden from the measurement image, but displayed as a "RUO spot".

Configuration


Calibration string:

Allows to override the calibration information from the cartridge's QR-Code. Leave empty to use the calibration information from the QR-Code on the cartridge (the usual behaviour).

Research spots:

Show research spots

If enabled, "research use only" spots will be displayed on measurement details and arrays. This does not affect the export behaviour or the reports of measurements.

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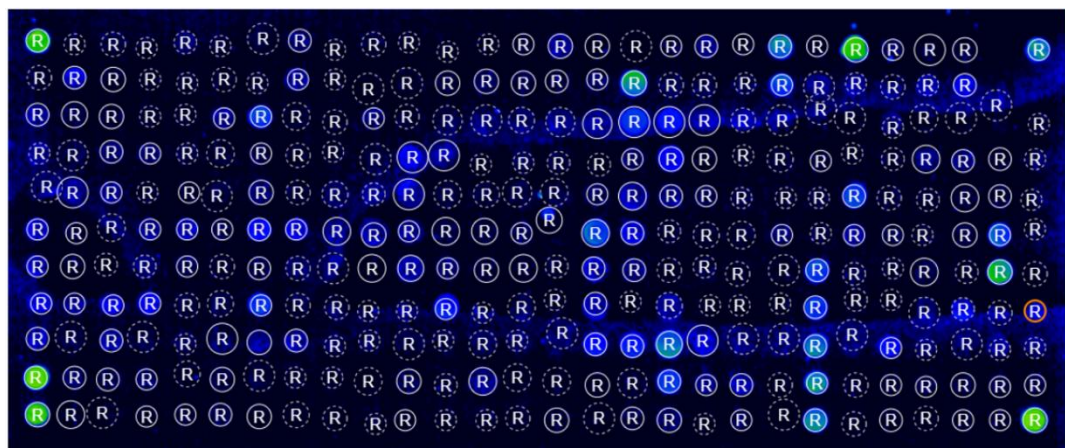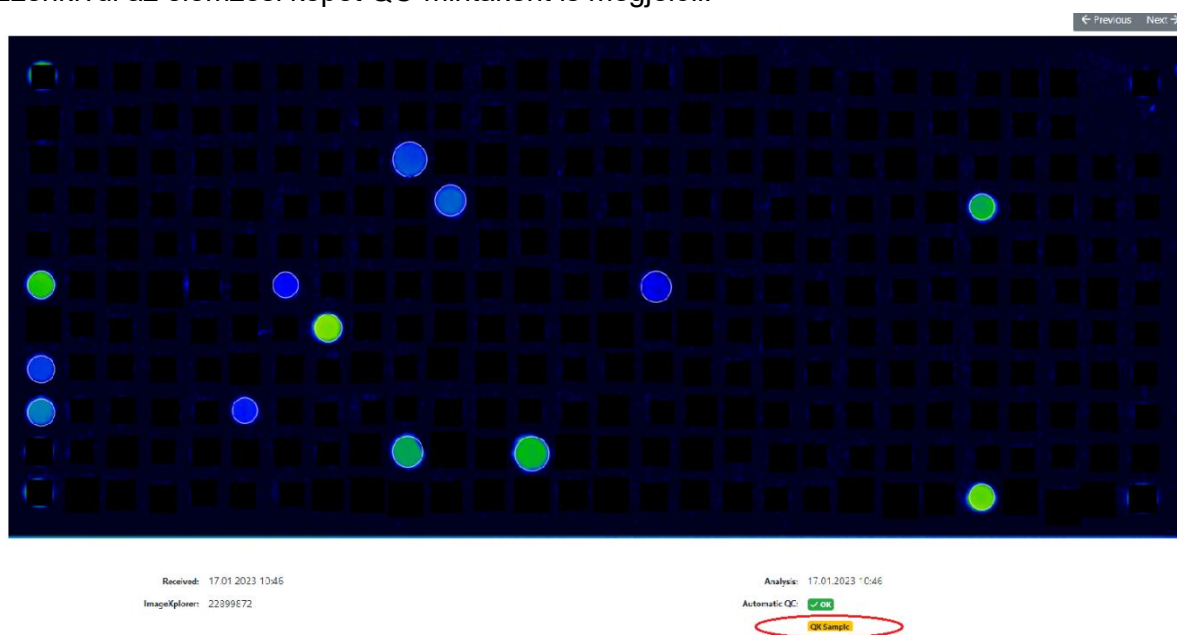
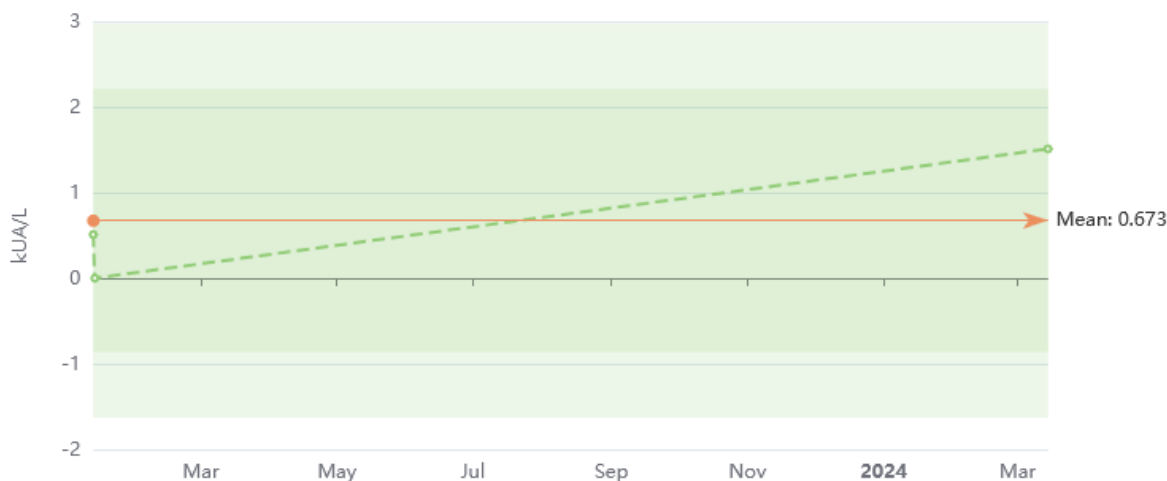
New Measurements

Search by barcode, patient, sample, tag, analysis date, a

Selected measurements ▾

- Download CSV format
- Download Excel format
- Approve
- Update interpretation
- Re-Process
- » Move
- Delete**

Tags	Product type
	FOX
	FOX
	ALEX2
	FOX
	ALEX2
	ALEX1

Delete measurements

You are about to delete the following measurements. Please note that this action is irreversible and cannot be undone.

Barcode	Analysis Date
02AAB160	16.12.2019 13:41

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Plant Food Food	M	Jug r 3	0.40
Plant Food Food	M	Jug r 1	0.38

Plant Food	M	Jug r 3	0.40
Food			
Plant Food	M	Jug r 1	0.38
Food			

Delete request scheduled now. [Cancel delete request](#)

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Sign in. Can be

Background jobs
View details about running background jobs and manually schedule runs.

[View dashboard](#)

Measurement delete / unapproval log & requests
Show logging and pending requests for measurement delete and unapproval.

[Show log & requests](#)

Measurement unapproval / delete requests

Open requests [Log & Finished requests](#)

Requested UTC	Requesting User (Email)	QR-Code	Request Type	Action
11.10.2023 13:05:28	Test Test (Test@test.com)	02AYQ145	Unapprove	Accept Reject
11.10.2023 12:59:27	Test Test (Test@test.com)	02AYQ14C	Delete	Accept Reject

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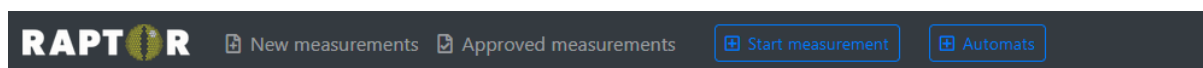
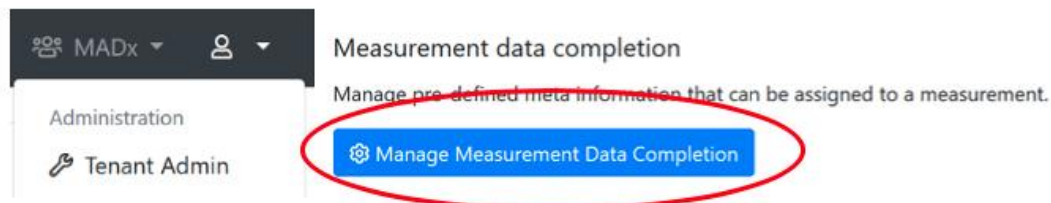
Add new user on Result Sharing Portal

General configuration

First name:

Last name:

E-Mail address:

Status:

Enabled (user can log in) 

Username:

Password:

Send an invitation e-mail to user automatically 

☒ Automatically send an email summary in case new measurement(s) have been shared
This email is sent once per day as a summary, just in case a new measurement was shared.

Doctor Info Text:


This information will be printed on a patients report.

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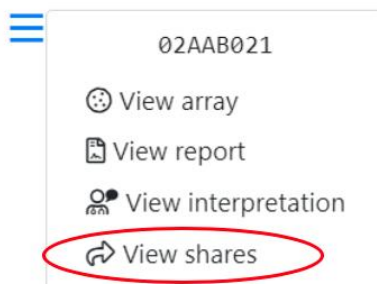≡ Share

You are about to share a measurement with a user of the User Portal. Please note that the doctor will have access to the results of the measurement. Once you share a measurement you'll be provided with a link to access the shared information via the encrypted User F

Barcode: 02AAB021

Share with: X Dr. (testuser2@macroarraydx.com)

If you can't find a specific user of the User Portal, go to the [tenant administration](#) to configure the list of available users of the User Portal.

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RAPTOR New measurements Approved measurements Start measurement with ImageXplorer Automats MADx Testtenant

≡ Share


Share this measurement

User	Link	Shared	
RP Bak	https://resultshare.raptor-server.com/s/07PXEpC7ZNyhm4PdRfmct4NR5pUtWq5FJcqOY7Ry6oEx08USj	7 months ago	 Revoke

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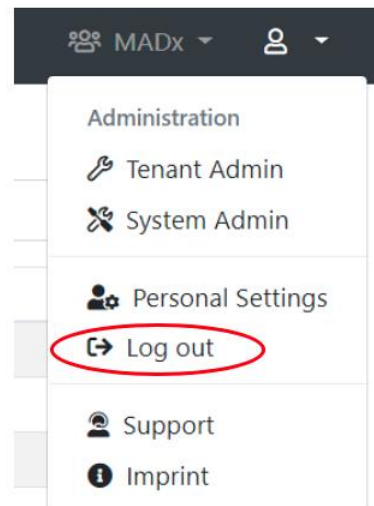Measurements

These are the measurements that are available for you. You may view the details of each measurement, or download the report PDF file.

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Anamneses

You can store anamneses for your patients, to assign them later to a measurement once it is available. An anamnesis then can be used to create a textual interpretation of the measurement.

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Edit anamnesis

Edit the details of an anamnesis here. This information is especially useful to know which anamneses is for what patient.

QR-Code E.g. 02ABC123	Sample Code E.g. 123456	Patient Code E.g. 987654	Patient Name E.g. Maria Musterfrau
Language English (en) ▼	Region Temperate ▼		
Notes E.g. Questionnaire filled by patient via our tablet.			
Save Cancel			

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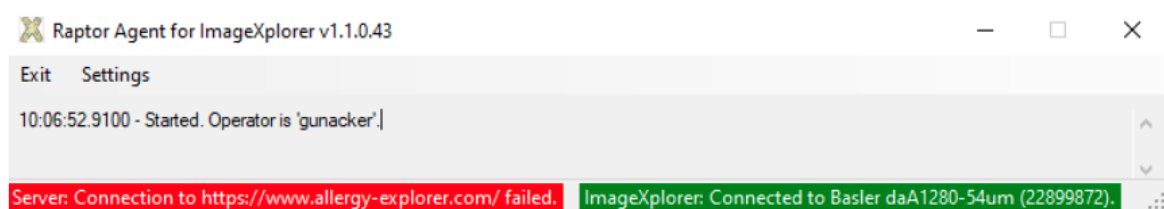
Settings ✕

Please enter the settings for your ImageXplorer. You can find your ImageXplorer key when you log in to Raptor. Please contact support if you're not sure what the proper settings are for your ImageXplorer.

Raptor URL:

ImageXplorer Key: 

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RAPTOR  [New measurements](#) [Approved measurements](#) [Start measurement with ImageXplorer](#)

The ImageXplorer IX training cannot be used, because it requires a ConfigXplorer scan.

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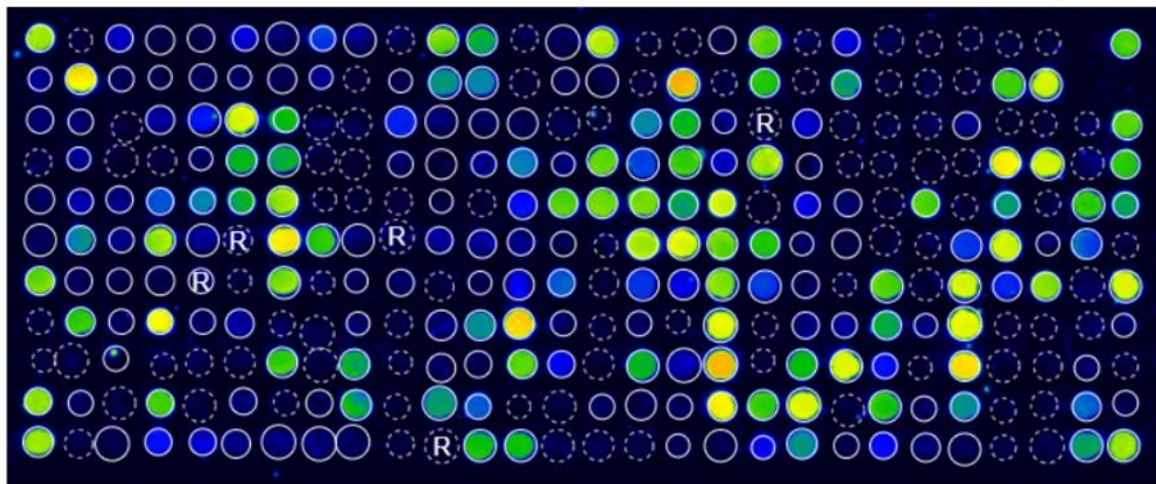
Analysis: 06.04.2023 13:37

ImageXplorer: 40287249

Automatic QC: ✗ FailedAnalysis Error: Insufficient data transfer - check connection.

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© Copyright by MacroArray Diagnostics

MacroArray Diagnostics (MADx)

Lemböckgasse 59, Top 4

1230 Vienna, Austria

+43 (0)1 865 2573

www.madx.com

Version number: 22-IFU-01-HU-10 (EN-27)

Date of Issue: 2026-08