# RAPTOR

SERVER ANALYSIS SOFTWARE

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### TARTALOMJEGYZÉK

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Welcome to RAPTOR SERVER Analysis Software, the analysis software for the [ALEX multiplex test](#).  
Please log in below to continue.

| Log in with Azure AD / Office 365                                                                                                                                                                                                                                                                       | Log in with Username & password                                                                                                                                                                                                                                                                                                                                                                                                      |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <p>If your organization integrates the software into its Microsoft Azure/Office 365 user directory (single-sign-on) click the button below. You will be redirected to the Microsoft login screen.</p> <p> Log in</p> | <p>Please enter your username and password.</p> <p><input type="text" value="Your username"/></p> <p><input type="password" value="Your password"/> </p> <p> Log in</p> <p>Forgot your password? No problem, you can <a href="#">request a new one</a>.</p> |

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Log in with Username & password

Please enter your username and password.

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**Contract approval required**

Please review the following documents. You are required to confirm that you have read and agree with each document in order to use RAPTOR SERVER:

☐ **Bitte lesen** (last updated 23.11.2023 18:00)

☐ **RAPTOR SERVER IFU** (last updated 04.12.2023 11:05)

☐ **Terms of Service** (last updated 22.11.2023 01:00)

Confirm

By clicking the button, I confirm that I have reviewed these documents and agree.

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| New Measurements                                                     |                       |          |             |                  | Last 7 days             |                        |
|----------------------------------------------------------------------|-----------------------|----------|-------------|------------------|-------------------------|------------------------|
| QC                                                                   | Product type / Source | QR-Code  | Sample Code | Analysis Date    |                         |                        |
| Failed                                                               | Research MAX 45k      | 99AFG0F7 | NK FOX      | 30/03/2023 20:50 | <a href="#">Details</a> | <a href="#">Report</a> |
| Failed                                                               | Research MAX 45k      | 99AFG0F8 | N.G.        | 30/03/2023 20:50 | <a href="#">Details</a> | <a href="#">Report</a> |
| Failed                                                               | Research MAX 45k      | 99AFG0F9 | E18735      | 30/03/2023 20:50 | <a href="#">Details</a> | <a href="#">Report</a> |
| Failed                                                               | Research MAX 45k      | 99AFG0FA | 20220325w56 | 30/03/2023 20:50 | <a href="#">Details</a> | <a href="#">Report</a> |
| Failed                                                               | Research MAX 45k      | 99AFG0EE | HM100038    | 30/03/2023 20:50 | <a href="#">Details</a> | <a href="#">Report</a> |
| Failed                                                               | Research MAX 45k      | 99AFG04D | 786118      | 30/03/2023 20:37 | <a href="#">Details</a> | <a href="#">Report</a> |
| Failed                                                               | Research MAX 45k      | 99AFG016 | HM100225    | 30/03/2023 20:37 | <a href="#">Details</a> | <a href="#">Report</a> |
| Failed                                                               | Research MAX 45k      | 99AFG0CA | E17685_M    | 30/03/2023 20:37 | <a href="#">Details</a> | <a href="#">Report</a> |
| Failed                                                               | Research MAX 45k      | 99AFG0CC | E17685_M    | 30/03/2023 20:37 | <a href="#">Details</a> | <a href="#">Report</a> |
| Failed                                                               | Research MAX 45k      | 99AFG015 | E17685_F    | 30/03/2023 20:37 | <a href="#">Details</a> | <a href="#">Report</a> |
| In total there are 270 new measurements in the selected time period. |                       |          |             |                  |                         |                        |
| <a href="#">View all new measurements</a>                            |                       |          |             |                  |                         |                        |


| Approved Measurements                          |                       |               |                  |                         |                        |  |
|------------------------------------------------|-----------------------|---------------|------------------|-------------------------|------------------------|--|
| QR-Code                                        | Product type / Source | Sample Code   | Analysis Date    |                         |                        |  |
| 02AVJ227                                       | ALEX2                 | Anna LMacey   | 28/03/2023 21:17 | <a href="#">Details</a> | <a href="#">Report</a> |  |
| 02ARA379                                       | ALEX2 IX              | SM_28032023   | 28/03/2023 14:29 | <a href="#">Details</a> | <a href="#">Report</a> |  |
| 02ARA393                                       | ALEX2 IX              | RG_28032023   | 28/03/2023 14:26 | <a href="#">Details</a> | <a href="#">Report</a> |  |
| 02AVA31F                                       | ALEX2 MAX 45k         | PR_Mar2023_11 | 21/03/2023 15:05 | <a href="#">Details</a> | <a href="#">Report</a> |  |
| 02AVI3B9                                       | ALEX2 MAX 45k         | KF            | 17/03/2023 14:46 | <a href="#">Details</a> | <a href="#">Report</a> |  |
| 02AYT337                                       | ALEX2 MAX 45k         | 2B46B46+CCD   | 14/03/2023 19:06 | <a href="#">Details</a> | <a href="#">Report</a> |  |
| 02ATL349                                       | ALEX2 IX              |               | 22/02/2023 14:21 | <a href="#">Details</a> | <a href="#">Report</a> |  |
| 02AVE383                                       | ALEX2 MAX 45k         | AO            | 03/02/2023 09:41 | <a href="#">Details</a> | <a href="#">Report</a> |  |
| 02AXI29A                                       | ALEX2 MAX 1k          | OW            | 19/01/2023 19:12 | <a href="#">Details</a> | <a href="#">Report</a> |  |
| 02AUY3AE                                       | ALEX2 MAX 45k         | BSZ_Okt2022_8 | 18/10/2022 18:12 | <a href="#">Details</a> | <a href="#">Report</a> |  |
| <a href="#">View all approved measurements</a> |                       |               |                  |                         |                        |  |

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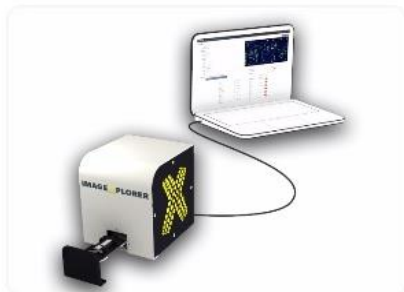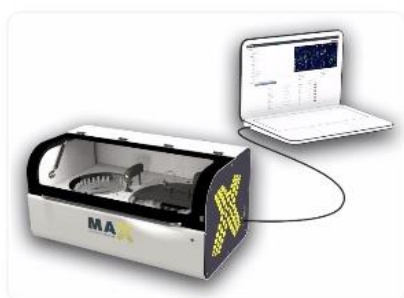
### 9K-Test-Instrument

Not initialized

[Initialize](#) [Shutdown](#) [Monthly Maintenance](#)

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**E-Mail address:**

**Language:** English (en)

The language is used to determine the display language.

**Culture:** English (United States) (en)

The culture is used to format numbers and dates.

**Timezone:** (UTC+01:00) Amsterdam, Berlin, Bern, Rome, Stockholm, Vienna

Configure your timezone to see date and times in your local time instead of UTC.

### Multi-Factor Authentication

[Manage MFA](#)

Use multi-factor authentication (MFA) to significantly increase security. It is highly recommended to enable MFA; this adds an additional layer of protection against phishing attacks or weak passwords.


### Configuration

**Research spots:** Show research spots

If enabled, "research use only" arrays. This does not affect the

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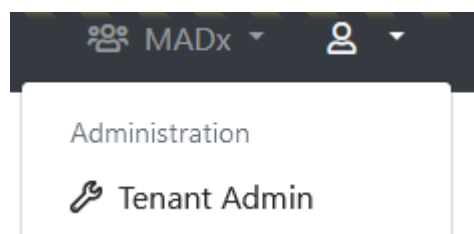### General configuration

**\* Name:** MADx

**Default timezone:** (UTC+01:00) Amsterdam, Berlin, Bern, R...  
Configure a timezone to see date and times in (instead of UTC). Please note that users can choose to use a different timezone.

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**Display MOD Excluded Spots As Research Use Only:** False  
If this setting is enabled, spots excluded by an MOD will not be hidden from the measurement image, but displayed as a "RUO spot".

**MFA:** ☐ Force MFA  
If this setting is enabled, users are forced to set up an MFA on sign-in, in case they haven't yet.

### MAX device default settings

**\* Default Minimum RunId:** 10  
When a new MAX device is created, this value will be used for its **Minimum Run Id**

### Tenant-specific batch releases

Restrict this tenant to only allow the following batches to be scanned or to restrict certain batches. Enter one per line.

Use:

☒ Whitelist  
☐ Blacklist

**Whitelisted Batches:**

Enter one batch code including product identifier per line, for example "02AAA"

**Blacklisted Batches:**

Enter one batch code including product identifier per line, for example "02AAA"

## Company Details

**\* Legal Entity Name:** MacroArray Diagnostics

**\* Street and Number:** Street 4

**\* Zip Code:** 1230

**\* City:** New York

**\* Country:** United States of America

**\* Phone:** +4318652573

**\* E-Mail:** office@macroarraydx.com

**\* Legal Representative:** Dr. Christian Harwanegg

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You requested a new measurement from ImageXplorer "IX Trainingsraum". Please wait until the measurement is completed ...

Invalid batch detected! The measurement can't be processed because its batch is blocked.

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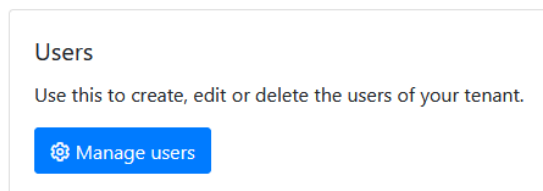
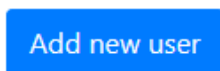
### Login Credentials

**Type:**

Please select if the user should be validated against Azure AD (Office 365), or if you want to specify username and password yourself. It is recommended to use Azure AD whenever possible. This allows for increased security by supporting two-factor authentication. Also, it is more convenient for the user, because he or she does not have to remember a new username and password (single-sign-on).

**Username:**

It is a good idea to use the users e-mail-address as his username, because the username needs to be unique across the entire system.

**Password:**

### Access Rights

#### Preconfigured groups

Use these pre-configured groups to automatically assign a set of access rights to this user. You can still add or remove specific access rights.

[Reset](#) [Lab User](#) [Super Lab User](#) [Support/Distributor User](#) [System Administrator](#)

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## Tenants

This user is connected to the following tenants and has the following access rights:

| <input type="checkbox"/> Tenants             | Can share measurements | Can edit measurements | Can edit personal data | Can approve measurements | Can manage basic tenant configuration | Can manage API keys | Can manage ImageExplorers | Can view agent connections | Can manage users | Can manage users of Result Sharing Portal | Can manage report templates | Can manage measurement data completion | Can unapprove tenant specific batch releases | Can use a custom calibration string | Can execute MAX d | Can skip t |
|----------------------------------------------|------------------------|-----------------------|------------------------|--------------------------|---------------------------------------|---------------------|---------------------------|----------------------------|------------------|-------------------------------------------|-----------------------------|----------------------------------------|----------------------------------------------|-------------------------------------|-------------------|------------|
| <input type="checkbox"/> Another MADx Tenant | ✓                      | ✓                     | ✓                      | ✓                        | ✓                                     | ✓                   | ✓                         | ✓                          | ✓                | ✓                                         | ✓                           | ✓                                      | ✓                                            | ✓                                   | ✓                 | ✓          |
| <input type="checkbox"/> MADx Tenant         | ✓                      | ✓                     | ✓                      | ✓                        | ✓                                     | ✓                   | ✓                         | ✓                          | ✓                | ✓                                         | ✓                           | ✓                                      | ✓                                            | ✓                                   | ✓                 | ✓          |

Additionally this user can be added to the following tenants by checking the tenant name and hitting save afterwards. User rights w applied from current open tenant user.

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Additionally this user can be added to the following tenants by checking the tenant name and hitting save afterwards. User right applied from current open tenant user.

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## Tenants

This user is connected to the following tenants and has the following access rights:

[Edit selected tenants rights](#)

| <input type="checkbox"/>            | Tenants             | Can share measurements              | Can edit measurements               | Can edit personal data              | Can approve measurements            | Can manage basic tenant configuration | Can manage API keys                 | Can manage ImageXplorers            |
|-------------------------------------|---------------------|-------------------------------------|-------------------------------------|-------------------------------------|-------------------------------------|---------------------------------------|-------------------------------------|-------------------------------------|
| <input checked="" type="checkbox"/> | Another MADx Tenant | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/>   | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> |

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Add new ImageXplorer

General information

|                            |                                                               |
|----------------------------|---------------------------------------------------------------|
| Name:                      | <input type="text"/>                                          |
| Key:                       | Cxltclg8okSwQ6CqJ4EZLWGnw2PdJ3                                |
|                            | <small>This key is required when installing an Agent.</small> |
| ImageXplorer Serialnumber: | <input type="text"/>                                          |
|                            | <small>Serialnumber of connected ImageXplorer.</small>        |
| Notes:                     | <div></div>                                                   |
| Internal Notes:            | <div></div>                                                   |
|                            | <small>For internal use only.</small>                         |
| Disabled:                  | Not disabled                                                  |

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
Details

Barcode: 02AE3488  
Sample Code: 12  
Requisition Code:  
Sampling Date: 01.10.2021  
Assay Date: 04.10.2021  
Patient Code: TE  
Patient Name: Test  
Date of Birth: 17.05.2010  
Info: Info  
Doctor Info: Dr. A  
Choose a doctor... Insert doctor info  
Tags: Tag Info  
MOD: < No MOD >  
Report Template: Version 2 (ALEX2 only) - Deutsch (de)  


Report Save Save and approve

← Previous Next →

Received: 04.10.2021 16:23  
ImageXplorer: 22124137  
Analysis: 04.10.2021 16:23  
Automatic QC: ✔ OK

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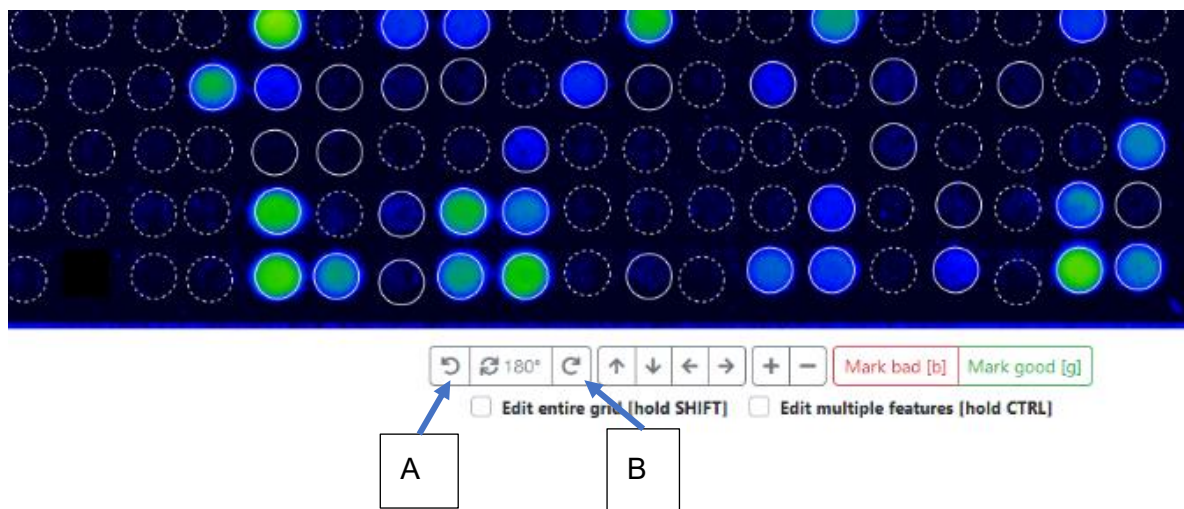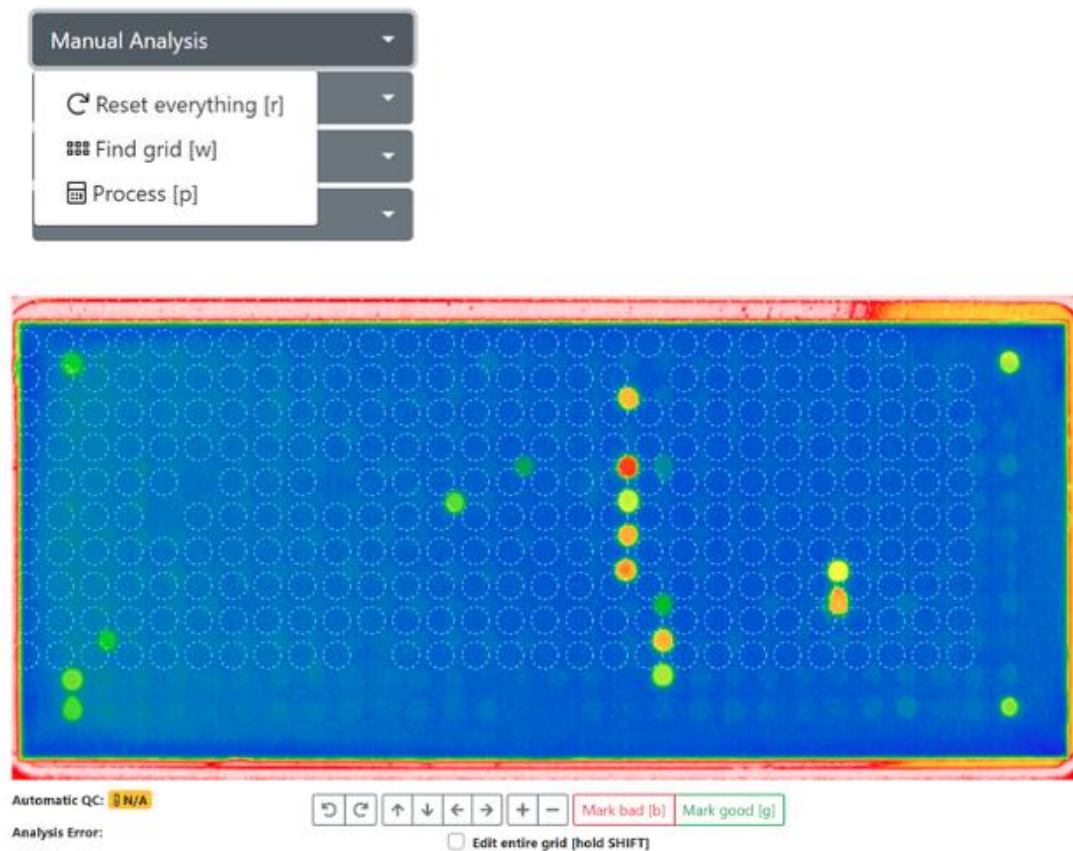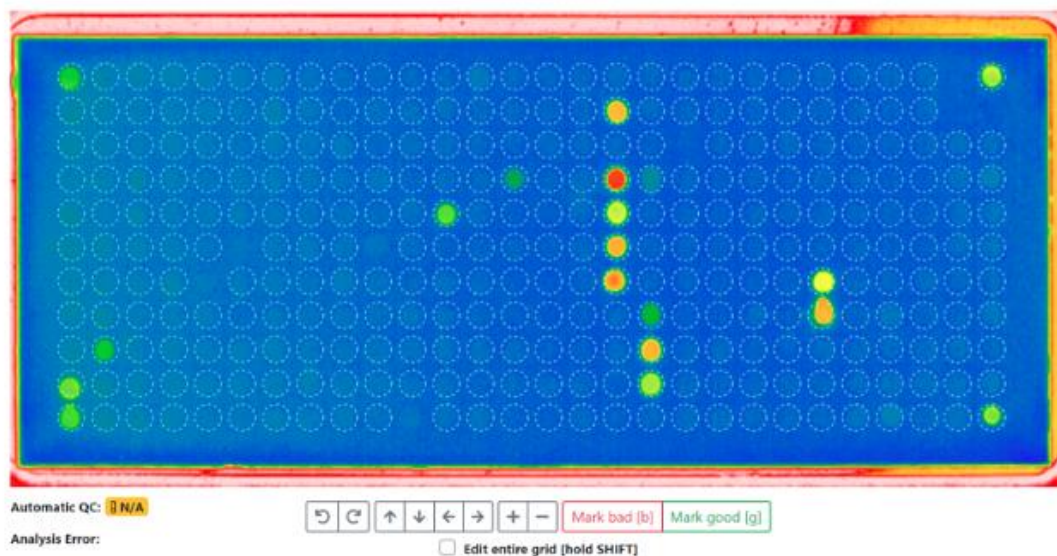
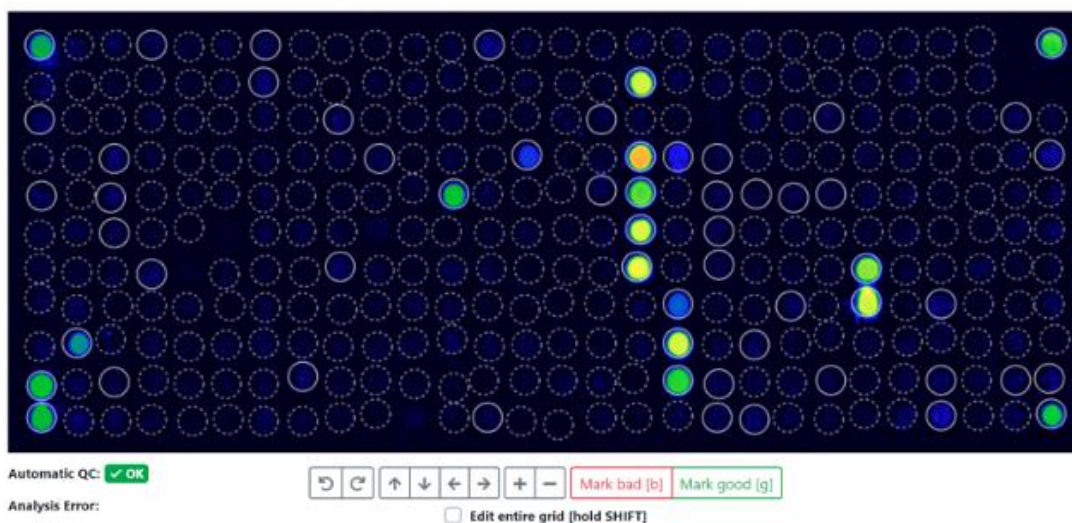
Manual Analysis
Image
Grid
More
Show Analysis & QC Logs

Reset everything and reprocess

This discards all modifications to the array, and re-runs the complete analysis process from the original image data. Please note that any manual changes will be irreversible lost.

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#### Analysis Log

```
'PreProcessedPixels' took 128ms.
Using definition '02_ALEXII_IVD_v1.6.0'.
'ApplyArrayDefinition' took 74ms.
Parsing calibration string ...
'ApplyDefinitions' took 75ms.
'FindGrid' took 1ms.
'AlignFeatures' took 862ms.
'CropImage' took 1ms.
'QuantifyFeatures' took 231ms.
'SubtractBackground' took 681ms.
Parsing calibration string ...
Using calibrators from QR-Code.
Total Igt (tIgt) normalized by calculating 0.330999755859375E-07x²+0.05697999954223633x+-38.69269921875
(x=8515), which results in 515.
'ReQuantify' took 57ms.
```

#### QC Log

These are the values used for QC (ALEX v2): Curve Correction=0.9997370867894242, tIgt=8515, GD1=28985 (SNR=105.38), GD2=29369 (SNR=185.28), GD3=28569 (SNR=147.53), GD4=28385 (SNR=128.1). Reference values are as follows: Guide Dots must have a value above 5000; each Guide Dot spot must exist and be unique; Guide Dot spot must not be set to bad. Average background must be between 4000 and 16000. QC (ALEX v2) succeeded.

#### Background Log

Calculated background statistics for 164 pixels (from a total of 182, skipped smallest 9, and highest 9 pixels). The average is 9,419, the minimum 7973 and the maximum 11792 (minimum and maximum is calculated including skipped pixels). This results in an average difference of 628, which is an SD of 49,17 and a CV of 0,52%.

For comparison, there is an overall total of 282944 pixels. The overall average is 14,871, the minimum 6,012 and the maximum 64,965.

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
| Tissue & Source      |  | Feature | kUA/L  |
|----------------------|--|---------|--------|
| Pollen<br>Inhalation |  | E Aca m | ≤ 0.10 |
| Mites                |  | E Aca s | 7.11   |

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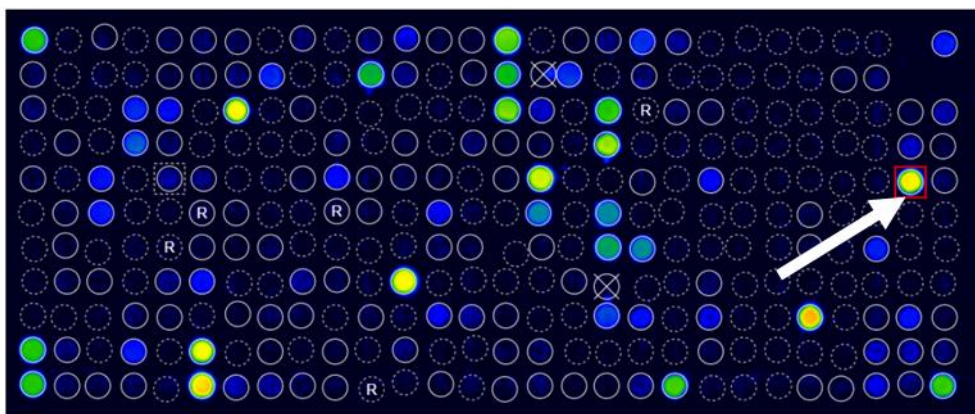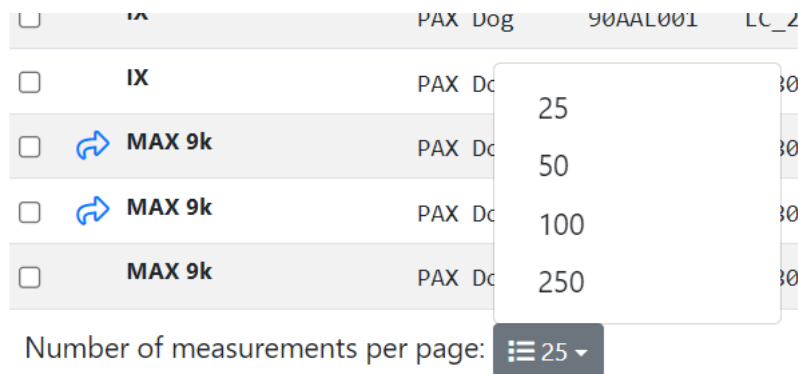
## New Measurements

Search by QR-Code, patient, sample, tag, analysis date, approval date ...

Selected measurements

| <input type="checkbox"/>            | Source | QC   | Tags | Product type |
|-------------------------------------|--------|------|------|--------------|
| <input type="checkbox"/>            | IX     | ✓ OK |      | ALEX2        |
| <input checked="" type="checkbox"/> | IX     | ✓ OK |      | ALEX2        |
| <input checked="" type="checkbox"/> | IX     | ✓ OK |      | ALEX2        |
| <input checked="" type="checkbox"/> | IX     | ✓ OK |      | ALEX2        |

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Delete

## Re-process measurements

You are about to re-process the following measurements. undone, the existing analysis result will be irreversibly over

| QR-Code  | Sample Code |
|----------|-------------|
| 02APK3C8 | Dummy8      |
| 02APK258 | Dummy1      |
| 02APK07F | Dummy13     |

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New Measurements

Search by QR-Code, patient, sample, tag, analysis date, approval date ...

|                          | Source | QC   | Tags | Product type | QR-Code  | Sample Code |
|--------------------------|--------|------|------|--------------|----------|-------------|
| <input type="checkbox"/> | IX     | ✓ OK |      | ALEX2        | 02AGT957 |             |
| <input type="checkbox"/> | IV     | ✓ OK |      |              |          |             |

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## New Measurements

Search by QR-Code, patient, sample, tag, analysis date, approval date ...

|                          | Source                              | QC   | Tags       | Product type |
|--------------------------|-------------------------------------|------|------------|--------------|
| <input type="checkbox"/> | Recaptured                          | ✓ OK | recaptured | ALEX1        |
| <input type="checkbox"/> | ⚠ Recaptured                        | ✓ OK | recaptured | ALEX1        |
| <input type="checkbox"/> | QR-Code was scanned more than once. |      | recaptured | ALEX1        |

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**RAPTOR** New measurements Approved measur

New Measurements

Search by QR-Code, patient, sample, tag, analysis date, approval date ...

Selected measurements ▾

|                                     | Source | QC   | Tags | Product type | QR-Code  |
|-------------------------------------|--------|------|------|--------------|----------|
| <input checked="" type="checkbox"/> | IX     | ✓ OK |      | ALEX2        | 02AMC180 |
| <input checked="" type="checkbox"/> | IX     | ✓ OK |      | ALEX2        | 02AMC078 |
| <input checked="" type="checkbox"/> | IX     | ✓ OK |      | ALEX2        | 02AKV2D9 |

Selected measurements ▾

- Download CSV format
- Download Excel format
- Approve**
- Update interpretation
- Re-Process
- » Move
- Delete

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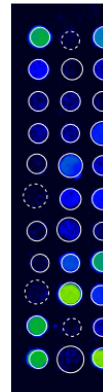
### Approve measurements

You are about to approve the following measurements.

There are measurements with one or more warnings. Please make sure to double-check these measurements before approval.

| QR-Code  | Sample Code | Analysis Date    |
|----------|-------------|------------------|
| 02AKV035 |             | 28.04.2022 11:14 |

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**RAPTOR**
New measurements
Approved measurements
Start measurement

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**RAPTOR**
New measurements
Approved measurements
Start measurement
Automats
MADx Testtenant

#### Exports

Certain types of data exports (for example, the export of many PDF reports at once) take a few minutes to prepare for download. Here download your requested exports. Please note that you can keep working, all exports are done in the background and you can download them when they are completed.

| Export requested   | # of measurements included | Export completed   |          |
|--------------------|----------------------------|--------------------|----------|
| 8/16/2021 11:07 AM | 4                          | 0 / 4              |          |
| 8/16/2021 11:06 AM | 4                          | 8/16/2021 11:06 AM |          |
| 8/16/2021 11:04 AM | 3                          | 8/16/2021 11:04 AM | Download |

Administration
Tenant Admin
Exports
Personal Settings
Log out
Support
Imprint
Download
Download

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| Export requested | # of measurements included |
|------------------|----------------------------|
| 06.07.2022 08:52 | 3                          |

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### Export multiple measurements

Please note that **Report PDF** and **ASTM** export is only a

|                                                    |                                              |                                       |
|----------------------------------------------------|----------------------------------------------|---------------------------------------|
| <b>QR-Code/Product/Batch</b>                       |                                              | <b>Sa</b>                             |
| <input type="text" value="02ABC"/>                 |                                              | <input type="text"/>                  |
| <b>Analysis Date (from)</b>                        | <b>Analysis Date (to)</b>                    |                                       |
| <input type="text" value="TT . MM . JJJJ"/>        | <input type="text" value="TT . MM . JJJJ"/>  |                                       |
| <b>Tag</b>                                         | <b>Status</b>                                |                                       |
| <input type="text"/>                               | <input type="text" value="All"/>             |                                       |
| <input type="checkbox"/> <b>Export CSV</b>         | <input type="checkbox"/> <b>Export Excel</b> |                                       |
| <input type="button" value="Search Measurements"/> |                                              | <input type="button" value="Cancel"/> |

## Export multiple measurements

Please note that **Report PDF** and **ASTM** export is only available for the current analysis.

QR-Code/Product/Batch

02aaa

Analysis Date (from)

TT . MM . JJJJ

Analysis Date (to)

TT . MM . JJJJ

Tag

Status

All

☐ Export CSV

☐ Export Excel

 Start export ...

← Back

66 Measurements found that match your criteria

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## Edit report template

Your changes have been saved successfully.

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PDF naming

**Setting Type:** Manual

The filename of the report PDF can be customized. Either built automatically based on your selection below or manually adjusted in the order of your needs.

**Name Pattern:** Report\_{PatientId}\_{MeasurementQrCode}

Available Elements: {PatientId} ; {PatientName} ; {SampleId} ; {MeasurementQrCode} ; {DateAndTime}

--

PDF naming

**Setting Type:** Automatic

The filename of the report PDF can be customized. Either built automatically based on your selection below or manually adjusted in the order of your needs.

**Include Patient ID:** ☒

**Include Patient Name:** ☒

**Include Sample ID:** ☐

**Include Measurement QR-Code:** ☐

**Include Date & Time:** ☒

Footer

**Footer text:**

Will be printed in the footer if provided and rendered as HTML.

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Configure report template "New report template"

The screenshot shows the 'Design' tab of a configuration interface. It has four sub-tabs: 'General', 'Design', 'Header', and 'Content'. The 'Design' tab is active. It contains three main settings:

- PDF format:** Two radio buttons, 'A4' (selected) and 'Letter'. A description to the right says: 'Specifies the paper size of the output PDF.'
- Font:** A dropdown menu showing 'Arial'. A description to the right says: 'The font family of the whole report texts can be adjusted.'
- Background image:** A text input field followed by a 'Browse' button. A description to the right says: 'If set, this image is used as background on the first page of the PDF.' Below the input field, it says 'Maximum file size 1 MB'.

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**Third level headline color:**

### Logos

**Logo (left):**  Browse  
Maximum file size 1 MB  
This logo will appear on the **top left corner** of the first page of the report.

**Text (center):**   
This text will appear either on the left, right, or center of the top of the first page (depending on your configured logos). Supports the placeholders `{{ PATIENT_PORTAL_QRCODE }}` and `{{ PATIENT_PORTAL_URL }}` to dynamically insert measurement information.

**Logo (right):**  Browse  
Maximum file size 1 MB  
This logo will appear on the **top right corner** of the first page of the report.

**Small logo:**  Browse  
Maximum file size 1 MB  
This logo is used on all detail pages. If not set, the default logo is used.

### Signatures

**Left:**  Browse  
Maximum file size 1 MB  
This signature or image will appear on the last page, left aligned.

**Center:**  Browse  
Maximum file size 1 MB  
This signature or image will appear on the last page, center aligned.

**Right:**  Browse  
Maximum file size 1 MB  
This signature or image will appear on the last page, right aligned.

### Indicators

**Custom intensity indicator:**   
A custom SVG for displaying allergen intensity can be pasted here. Will be styled according to specified colors.

**⚠** The SVG element must contain the attributes **width** and **height** with a minimum dimension of 16 for correct display in PDF.

**🌿** SVGs paths which should be colored must require the **colorize** class as well as a fill defined in the **defs** segment of the SVG.

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General Design Header **Content**

Header content:

| Left half                  |                                     | Right half           |                          |
|----------------------------|-------------------------------------|----------------------|--------------------------|
| Custom headline            |                                     | Custom headline      |                          |
| Patient Name               | <input type="checkbox"/>            | Doctor Information   | <input type="checkbox"/> |
| Patient Name               | <input checked="" type="checkbox"/> | None                 | <input type="checkbox"/> |
| Custom headline            |                                     |                      |                          |
| Date of Birth              | <input type="checkbox"/>            | None                 | <input type="checkbox"/> |
| Sample ID                  | <input type="checkbox"/>            | None                 | <input type="checkbox"/> |
| QR-Code                    | <input type="checkbox"/>            | Additional Info text | <input type="checkbox"/> |
| Analysis Date              | <input type="checkbox"/>            | None                 | <input type="checkbox"/> |
| Amount of tested allergens | <input type="checkbox"/>            | None                 | <input type="checkbox"/> |
| Test method                | <input type="checkbox"/>            | Quality Control text | <input type="checkbox"/> |

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| General                                                     | Design                              | Header                                                                                                                                                                    | Content |
|-------------------------------------------------------------|-------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------|
| <b>Allergens to include:</b>                                | All allergens                       | Specifies which kind of allergens should be included in the report.                                                                                                       |         |
| <b>Allergens sort order:</b>                                | Alphabetically by allergen          | Specifies the sort order of the allergens.                                                                                                                                |         |
| <b>Include Overview:</b>                                    | <input checked="" type="checkbox"/> | Whether the allergen overview should be printed or not.                                                                                                                   |         |
| <b>Include product type in overview:</b>                    | <input checked="" type="checkbox"/> | Whether the product type should be printed in the header of allergen overview.                                                                                            |         |
| <b>Include Total IgE:</b>                                   | On the front page                   | Whether the information of total IgE should be included in the report and where.                                                                                          |         |
| <b>Cross-reactive information:</b>                          | Include text                        | This section will be displayed underneath the allergen list itself.                                                                                                       |         |
| <b>Include Interpretation:</b>                              | Include manually                    | It can be specified if the Raven report interpretation information should be included in the report or not.<br>This section is printed underneath the allergens overview. |         |
| <b>i Additional information (below the list of results)</b> |                                     | Those informations are integrated into the report in an extra section underneath the allergen list.                                                                       |         |
| <b>Approval date:</b>                                       | <input type="checkbox"/>            |                                                                                                                                                                           |         |
| <b>Assay date:</b>                                          | <input type="checkbox"/>            |                                                                                                                                                                           |         |
| <b>Doctor info:</b>                                         | <input type="checkbox"/>            | The information about the doctor which accompanied the measurement/patient.                                                                                               |         |
| <b>Print date:</b>                                          | <input checked="" type="checkbox"/> | When the report was generated.                                                                                                                                            |         |
| <b>Report info:</b>                                         | <input type="checkbox"/>            | The „additional info” field belonging to the measurement.                                                                                                                 |         |
| <b>Sample date:</b>                                         | <input checked="" type="checkbox"/> |                                                                                                                                                                           |         |

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#### Configure report template "New report template - Duplicate"

| General                    | Design                          | Header | Content |
|----------------------------|---------------------------------|--------|---------|
| <b>Name:</b>               | New report template - Duplicate |        |         |
| <b>Support products:</b>   | ALEX <sup>2</sup>               |        |         |
| <b>Preferred language:</b> | Deutsch (de)                    |        |         |

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**General configuration**

**\* Name:** MADx

**Default timezone:** (UTC+01:00) Amsterdam, Berlin, Bern, Rome, Stockholm, Vienna  
Configure a timezone to see date and times in (instead of UTC). Please note that users can choose to use a different timezone.

**Default report template (ALEX 1):** Version 1 (ALEX1 & ALEX2)

**Default report template (ALEX 2):** Version 2 (ALEX2 only)

**Default report language:** English (United Kingdom) (en-gb)

**Default MOD:** < No MOD >

**Require ImageXplorer testarray scan:** Warn after 60 days, enforce after 90 days.

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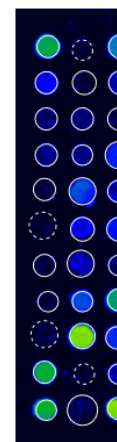
Template: Version 2 (ALEX2 only)

Language: English (en)

Usually this tenant's default template and language will be used when sharing or exporting measurements. You may change the template and language associated with this measurement by saving your current configuration.

You can download this report in PDF format once the measurement has been approved.

You have access to manage the report templates and settings via the tenant administration. Please note that any changes will affect all users.

**MADX**  
MACROARRAY DIAGNOSTICS

PATIENT ID  
PATIENT NAME  
DATE OF BIRTH  
SAMPLE ID  
ASSAY CODE  
02APH008  
ANALYZED ON  
11/04/2022  
TESTED ALLERGENS  
202  
TEST METHOD CODE  
ALEX2

REFERRING PHYSICIAN

NOTE  
The internal QC (Fluorescence check for QC) was within acceptance range.

Lab report: Summary on detectable sensitisations

POLLEN  
Grass Pollen

MICROORGANISMS  
Fungal Spores & Yeast

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**RAPTOR** New measurements Approved measurements Start m

**Approved Measurements**

Search by QR-Code, patient, sample, tag, analysis date ...

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**Tags:**

**MOD:**

**Report Template:** Version 2 (ALEX2 only) - Deutsch (de)

[New measurements](#)
[Approved measurements](#)
[Start measurement](#)
[Automats](#)

### Report

Template: 
Language:

[Save](#)
[Use defaults](#)

Usually this tenants default template and language will be used when sharing or exporting measurements. You may change the template and language associated with this measurement by saving your current configuration.

[Download PDF](#)

[Manage report template and settings](#)

You have access to manage the report templates and settings via the tenant administration. Please note that any changes will affect all users.

|                |                                |
|----------------|--------------------------------|
| Patient ID:    |                                |
| Patient Name:  |                                |
| Date of Birth: |                                |
| Sample ID:     |                                |
| Barcode:       | 02AAE                          |
| Sampled on:    |                                |
| Tested on:     | 12/18/                         |
| Printed on:    | 12/19/                         |
| Note:          | The internal QC (Plausibility) |

Lab report: Summary on de


|                  |               |
|------------------|---------------|
| Pollen           | Grass Pollen  |
|                  | Tree Pollen   |
|                  | Weed Pollen   |
| Mites            | House Dust M  |
| Microorganisms   | Fungal Spores |
|                  | Legumes       |
|                  | Grain         |
| Plant-Based Food | Spices        |
|                  | Fruits        |
|                  | Vegetables    |
|                  | Nuts & Seeds  |
|                  | Milk          |

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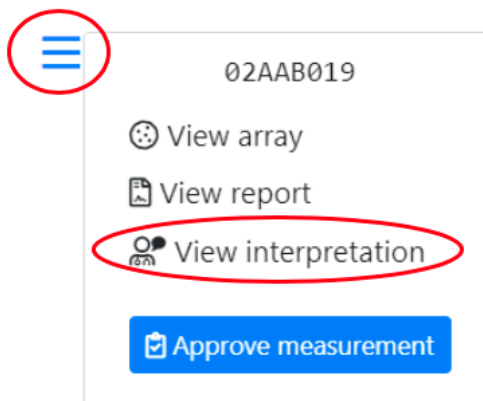**Rule set**

RAVEN Ruleset A for ALEX

A rule set defines the RAVEN or RAVEN<sup>2</sup> engine that is used to interpret results.

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AIT

Medication

### MADx data declaration for smartRaven

☒ I accept to enter clinical data for using the full interpretation function of smartRAVEN

### Patient has reported any symptoms

☐ Yes ☐ No

### Additional information about the patient

☐ Atopic Dermatitis

☐ Allergen-specific immunotherapy (AIT) - ongoing, completed or discontinued

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[Request new interpretation](#)

This measurement already has an interpretation stored. Usually, it should not be necessary to request those again, but it could be necessary when the patients information has changed.

#### Interpretation

##### Raven Interpretation Summary

##### Sample Information

The sample was tested on ALEX<sup>2</sup> Barcode 02AAB019, interpretation date 12/19/2019.

Of the tested 295 allergens, 25 were above the detection limit of 0.3 kU<sub>L</sub>/L. A sensitization can be an indicator of an allergy. For selected allergens, comments for interpretation are listed below.

##### Total IgE: 512 kU/L

The measured total IgE was 512 kU/L. Above a total IgE of 100 kU/L, allergy is likely.

##### Cross-Reactive Allergen Sensitization detected

Sensitization against allergen molecules which are markers of (broad) cross-reactivity between different allergen sources.

Detected cross-reactive allergen sensitizations:

- PR-10s: Aln g 1, Bet v 1, Cor a 1.0103, Fag s 1
- Polcalcins: Aln g 4, Phi p 7

##### PR-10 Proteins (PR10)

PR-10 inhalative: The major birch pollen allergen, Bet v 1, represents the prototype of all PR-10 allergens and is the primary sensitizer in regions with birch pollen exposure. The presence of PR-10 allergens in Fagales tree pollen explains IgE cross-reactivity between pollen from hazel, alder, beech, oak and hornbeam. PR-10 nutritive: PR-10 allergens in raw fruits, nuts, vegetable and legumes can induce oral allergy syndrome and sometimes severe allergic reactions in sensitized individuals. PR-10 allergens are not stable to heat and digestion.

##### Polcalcins (Ca++ binding proteins, PC)

Members of the PC family can cause inhalative symptoms. PC allergens can be found in pollen from grasses, trees and weeds. Inhalative symptoms manifest as allergic rhino-conjunctivitis and/or allergic asthma.

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The screenshot displays the 'Details' form in the MacroArray software. The form includes fields for Barcode (024E1852), Sample Code, Sampling Date (tt.mm.jjj), Assay Date (07.10.2020), Patient Code, Patient Name, Date of Birth (tt.mm.jjj), Info, Doctor Info, Tags, and MODs (Pollen allergens). A red oval highlights the MODs field. Below the form are buttons for Report, Save, and Save and approve. The bottom section shows the resulting spot array image, which is a grid of spots with blue and green fluorescence. Navigation buttons (Previous, Next) are visible above the image. At the bottom, there is a summary bar with the following information:

| Received:              | Analysis:                                             |
|------------------------|-------------------------------------------------------|
| 07.10.2020 14:44       | 07.10.2020 14:44                                      |
| ImageXplorer: 23156891 | Automatic QC: <span style="color: green;">✔ OK</span> |

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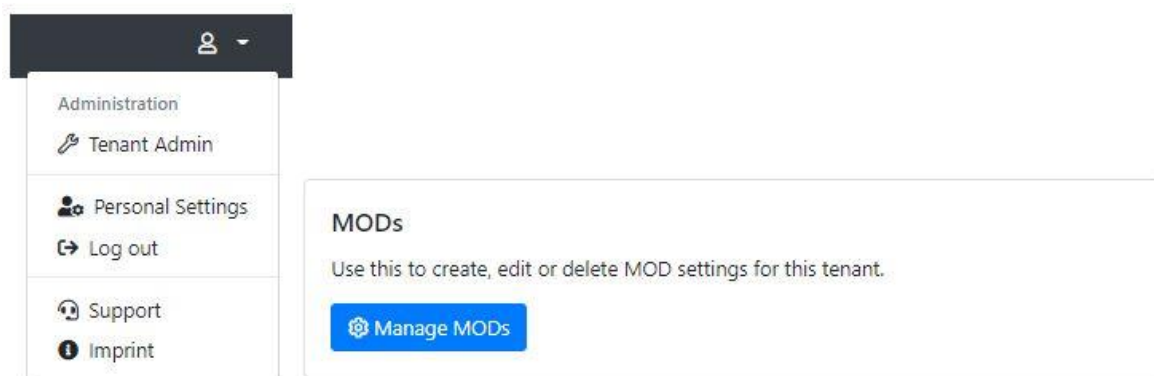
**Type:**

**Deactivated:**

You cannot delete an MOD once it is in use. But you can deactivate it, so that it's not available for new measurements anymore.

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**\* Name:**

**Default timezone:**

Configure a timezone to see date and times in (instead of UTC). Please note that users can choose to use a different timezone.

**Default report template (ALEX 1):**

**Default report template (ALEX 2):**

**Default report template (FOX):**

**Default report template (PAX):**

**Default report language:**

**Default ALEX2 MOD:**

**Default FOX MOD:**

**Default PAX MOD:**

**Force MOD:**

Requires a measurement to always have an MOD assigned. If this setting is enabled, a measurement cannot be approved without an MOD assigned.

**Hide Global MOD:**

If this setting is enabled, only the tenant-specific MOD will be available to select, and the global system-wide MOD will be hidden.

**Display MOD Excluded Spots As Research Use Only:**

If this setting is enabled, spots excluded by an MOD will not be hidden from the measurement image, but displayed as a "RUO spot".

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MOD:

Report Template:

**Force MOD:**   
 Requires a measurement to always have an MOD assigned. If this setting is enabled, a measurement cannot be approved without an MOD assigned.

**Hide Global MOD:**   
 If this setting is enabled, only the tenant-specific MOD will be available to select, and the global system-wide MOD will be hidden.

**Display MOD Excluded Spots As Research Use Only:**   
 If this setting is enabled, spots excluded by an MOD will not be hidden from the measurement image, but displayed as a "RUO spot".

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**Force MOD:** False  
Requires a measurement to always have an MOD assigned. If this setting is enabled, a measurement cannot be approved without an MOD assigned.


**Hide Global MOD:** False  
If this setting is enabled, only the tenant-specific MOD will be available to select, and the global system-wide MOD will be hidden.

**Display MOD Excluded Spots As Research Use Only:** True  
If this setting is enabled, spots excluded by an MOD will not be hidden from the measurement image, but displayed as a "RUO spot".

**Configuration**

**Calibration string:**   
Allows to override the calibration information from the cartridge's QR-Code. Leave empty to use the calibration information from the QR-Code on the cartridge (the usual behaviour).

**Research spots:** Show research spots  
If enabled, "research use only" spots will be displayed on measurement details and arrays. This does not affect the export behaviour or the reports of measurements.

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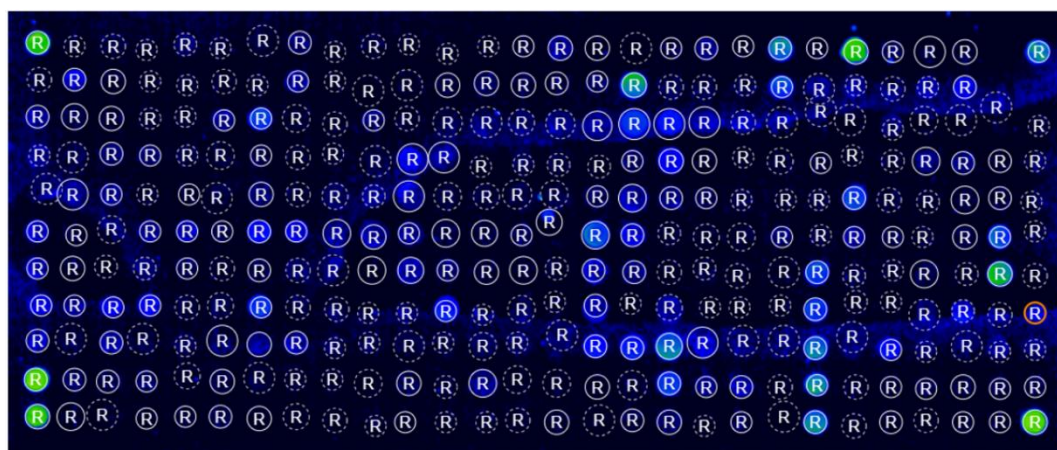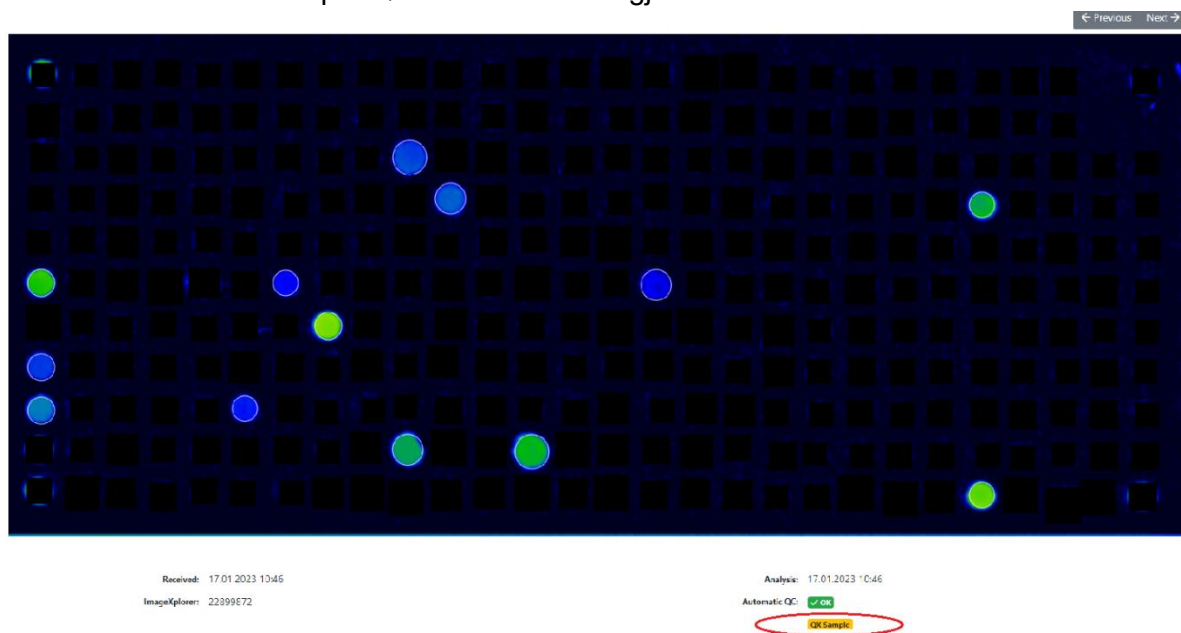
### New Measurements

Search by barcode, patient, sample, tag, analysis date, a

Selected measurements ▾

- Download CSV format
- Download Excel format
- Approve
- Update interpretation
- Re-Process
- » Move
- Delete**

| Tags | Product type |
|------|--------------|
|      | FOX          |
|      | FOX          |
|      | ALEX2        |
|      | FOX          |
|      | ALEX2        |
|      | ALEX1        |

### Delete measurements

You are about to delete the following measurements. Please note that this action is irreversible and cannot be undone.

| Barcode  | Analysis Date    |
|----------|------------------|
| 02AAB160 | 16.12.2019 13:41 |

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|                 |   |         |      |
|-----------------|---|---------|------|
| Plant Food Food | M | Jug r 3 | 0.40 |
| Plant Food Food | M | Jug r 1 | 0.38 |

|            |   |         |      |
|------------|---|---------|------|
| Plant Food | M | Jug r 3 | 0.40 |
| Food       |   |         |      |
| Plant Food | M | Jug r 1 | 0.38 |
| Food       |   |         |      |

Delete request scheduled now. [Cancel delete request](#)

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Sign in. Can be

**Background jobs**  
View details about running background jobs and manually schedule runs.

[View dashboard](#)

**Measurement delete / unapproval log & requests**  
Show logging and pending requests for measurement delete and unapproval.

[Show log & requests](#)

#### Measurement unapproval / delete requests

Open requests

[Log & Finished requests](#)

| Requested UTC       | Requesting User (Email)   | QR Code  | Request Type | Action                                        |
|---------------------|---------------------------|----------|--------------|-----------------------------------------------|
| 11.10.2023 13:05:28 | Test Test (Test@test.com) | 02AYQ145 | Unapprove    | <div><div>Accept</div><div>Reject</div></div> |
| 11.10.2023 12:59:27 | Test Test (Test@test.com) | 02AYQ14C | Delete       | <div><div>Accept</div><div>Reject</div></div> |

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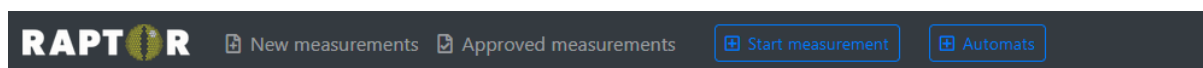
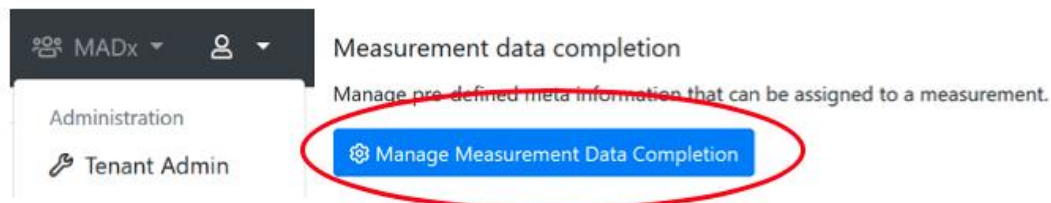
### Add new user on Result Sharing Portal

General configuration

First name:

Last name:

E-Mail address:

Status:

Enabled (user can log in) 

Username:

Password:

Send an invitation e-mail to user automatically 

☒ Automatically send an email summary in case new measurement(s) have been shared  
This email is sent once per day as a summary, just in case a new measurement was shared.

Doctor Info Text:


This information will be printed on a patients report.

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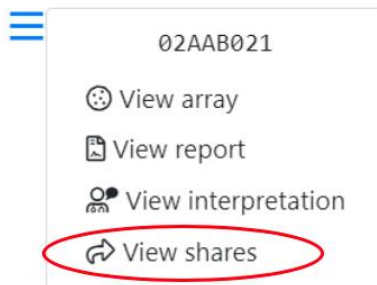## Share

You are about to share a measurement with a user of the User Portal. Please note that the doctor will have access to the results of the measurement. Once you share a measurement you'll be provided with a link to access the shared information via the encrypted User F

Barcode: 02AAB021

Share with: X Dr. (testuser2@macroarraydx.com)

If you can't find a specific user of the User Portal, go to the [tenant administration](#) to configure the list of available users of the User Portal.

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**RAPTOR** New measurements Approved measurements Start measurement with ImageXplorer Automats MADx Testtenant

Share


Share this measurement

| User   | Link                                                                                                                                                                                              | Shared       |        |
|--------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------|--------|
| RP Bak | <a href="https://resultshare.raptor-server.com/s/07PXEpC7ZNyhm4PdRfmct4NR5pUtWq5FJcqOY7Ry6oEx08USj">https://resultshare.raptor-server.com/s/07PXEpC7ZNyhm4PdRfmct4NR5pUtWq5FJcqOY7Ry6oEx08USj</a> | 7 months ago | Revoke |

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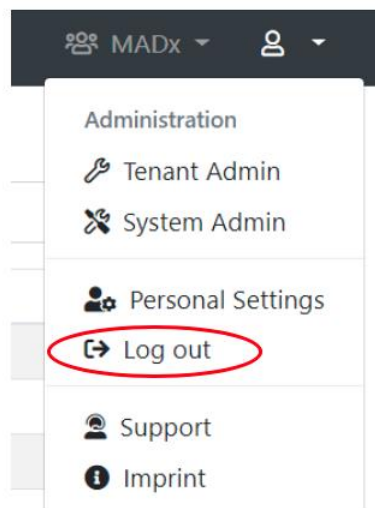 Result Sharing Portal English (en) ▼

### Welcome to the Result Sharing Portal

The Raptor-Server Result Sharing Portal allows users to view shared measurements. You will need a username and password to log in to view all measurements that have been shared with you. Also, a sharing link will allow to view specific shared measurements.

Please enter your username and your password:

Did you forget your password? No problem, you can [request a new one](#).

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Result Sharing Portal
English (en) ▾
Change password
Log out

### Measurements

These are the measurements that are available for you. You may view the details of each measurement, or download the report PDF file.

| Measurement              | Product | Patient | Sample   | Date       |
|--------------------------|---------|---------|----------|------------|
| <a href="#">01AAK239</a> | ALEX1   |         | 11852162 | 7/12/2019  |
| <a href="#">01AAK281</a> | ALEX1   |         | 04824821 | 12/18/2019 |
| <a href="#">01AAK891</a> | ALEX1   |         | 66560260 | 1/28/2020  |
| <a href="#">02AAA015</a> | ALEX2   |         |          | 2/21/2020  |

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## Measurement 02AYQ143

Product: ALEX2  
Patient:  
Date of Birth:  
Sample:  
Info:  
Doctor Info:  
Sampling Date:  
Assay Date: 10/12/2023  
Processed Date: 10/12/2023  
Approval Date: 10/12/2023

**Report Language:**

English (United Kingdom) ▼

You can request a new RAVEN interpretation that will be included in the report.

[Request new RAVEN interpretation](#)

You can download a PDF report that includes all information and results of this measurement.

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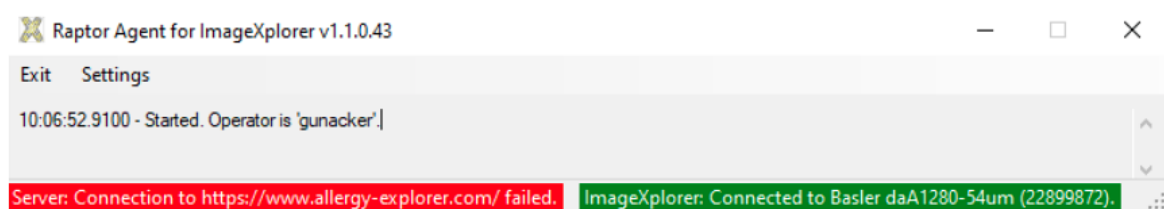
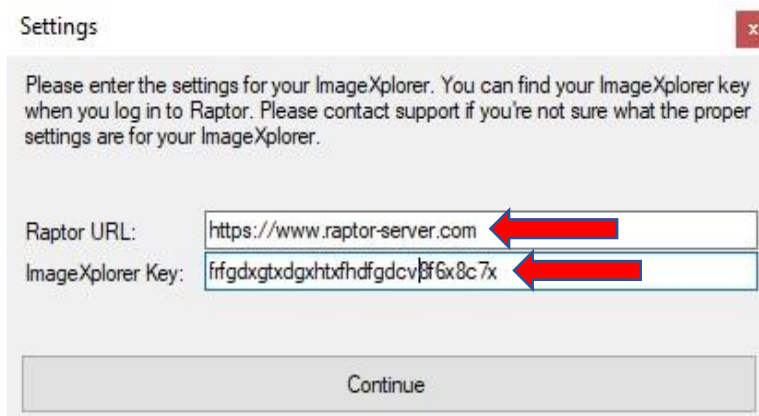
**Settings** ✕

Please enter the settings for your ImageXplorer. You can find your ImageXplorer key when you log in to Raptor. Please contact support if you're not sure what the proper settings are for your ImageXplorer.

Raptor URL:

ImageXplorer Key:  

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**RAPTOR**  [New measurements](#) [Approved measurements](#) [Start measurement with ImageXplorer](#)

The ImageXplorer IX training cannot be used, because it requires a ConfigXplorer scan.

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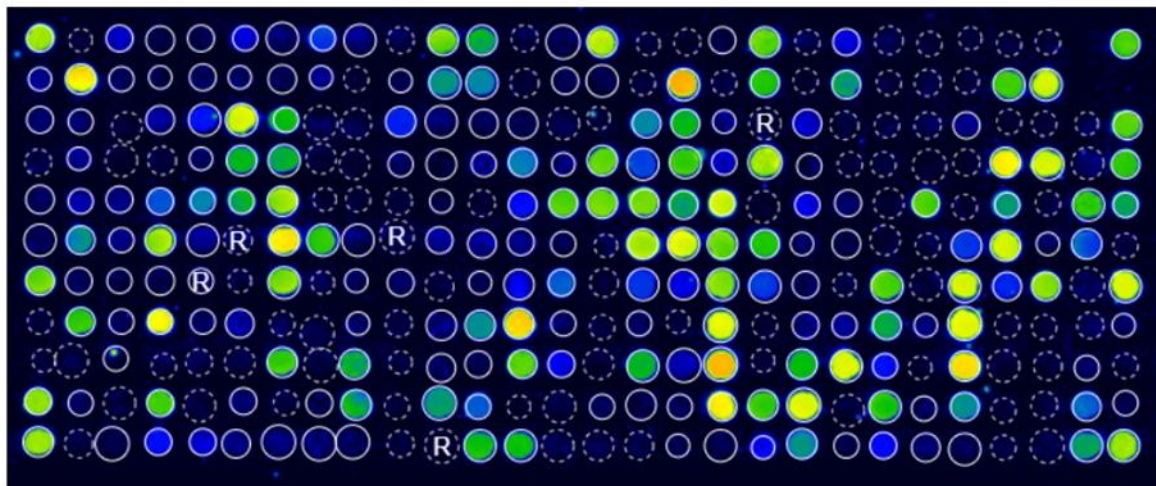
Analysis: 06.04.2023 13:37

ImageXplorer: 40287249

Automatic QC: ✖ FailedAnalysis Error: Insufficient data transfer - check connection.

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© Copyright by MacroArray Diagnostics

MacroArray Diagnostics (MADx)

Lemböckgasse 59, Top 4

1230 Vienna, Austria

+43 (0)1 865 2573

[www.madx.com](http://www.madx.com)

Version number: 22-IFU-01-HU-09 (EN-26)

Date of Issue: 2026-02