



BeGraft⁺

Now
BEVAR
indicated

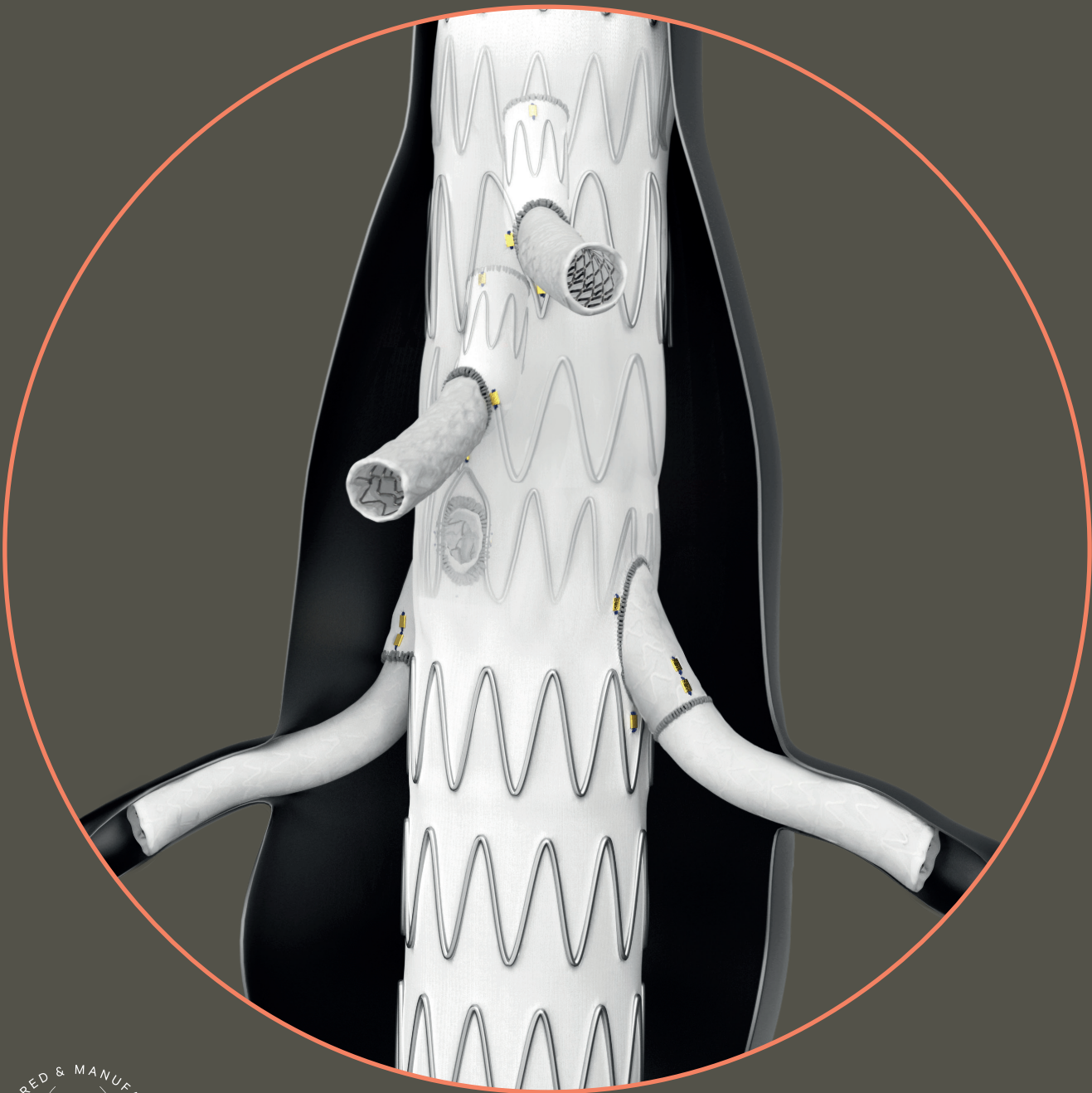

Built for BEVAR. Built to last.

Bridging distances in BEVAR interventions require heavy support due to movement and angulation of the target vessels. The BeGraft Plus is engineered for BEVAR procedures and is compatible with inner- and outer branches and with custom-made and off-the-shelf-solutions. It complements Bentley's broad portfolio of solutions for complex aortic bridging interventions that are dedicated to patient needs.

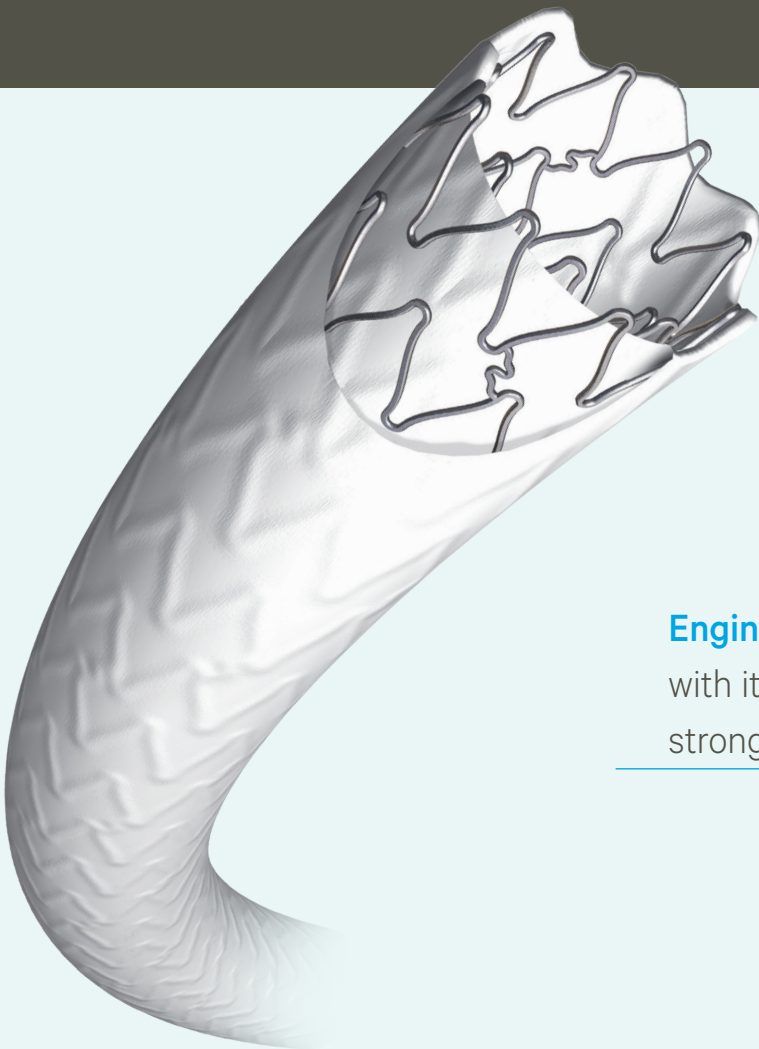

Bridging of longer distances in BEVAR

due to lengths up to 78 mm

Engineered to reduce the risk of endoleaks

with its unique stent-in-stent design that provides strong fixation and secure sealing


Excellent patency rates as bridging stent¹

as well as highest radial force and exceptional kink resistance² while still maintaining flexibility in lateral direction

¹The BeGraft Plus is based on the design of the BeGraft Peripheral Plus, see Trial CT11-01

²compared to other available covered stents, data on file

THE DOUBLE ONE

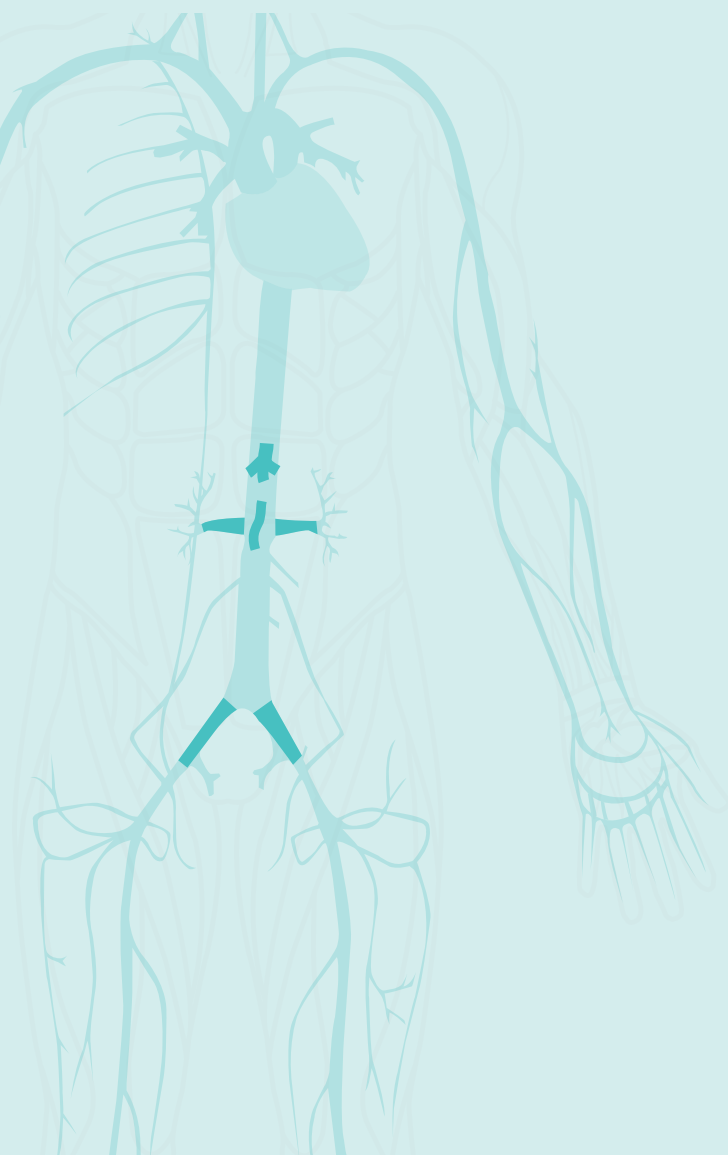

Product specifications

Stent Material	CoCr (L-605)
Graft Material	Two layers of micro-porous ePTFE tubing (203 ± 25µm)
Stent Graft Design	Double stent scaffold Open-cell design Patented ePTFE cover fixation
Delivery System	0.035" OTW
Catheter Shaft Length	120 cm
Expanded Stent Graft Diameter	Ø 5 and 6 mm ² Ø 7 and 8 mm ³ Ø 9 and 10 mm ⁴
Nominal Stent Graft Length	28,38,58,78 mm ² 27,37,57,77 mm ^{3,4}
Introducer Sheath Compatibility	7/8 F ⁵ for Ø 5 - 7 mm 8 F for Ø 8 mm 8/9/10 F ⁵ for Ø 9 and 10 mm
Shelf Life	3 years

²SV= Small Vessel Design, ³MV= Medium Vessel Design,

⁴LV= Large Vessel Design,

⁵ Depending on indication and access site

Three different stent designs ensure ideal stent performance.

Indications

The BeGraft Plus is indicated for permanent intraluminal placement in iliac and renal arteries for

- restoring and improving the patency impaired by arteriosclerosis or other stenoses that may be treated by stenting
- sealing aneurysms, acute perforations, acute ruptures, and fistulas.

Furthermore, the BeGraft Plus is indicated for use

- as bridging stent graft in patients requiring endovascular abdominal aortic aneurysm repair with a branched aortic endovascular graft to keep the visceral arteries perfused.

Ordering information

Expanded Stent Graft Diameter	Nominal Stent Graft Length	Introducer Sheath Size	Catalogue Number for Catheter Length 120 cm
5 mm	28 mm	7/8 F ⁵	BGP+2805_2
	38 mm		BGP+3805_2
	58 mm		BGP+5805_2
	78 mm		BGP+7805_2
6 mm	28 mm	7/8 F ⁵	BGP+2806_2
	38 mm		BGP+3806_2
	58 mm		BGP+5806_2
	78 mm		BGP+7806_2
7 mm	27 mm	7/8 F ⁵	BGP+2707_2
	37 mm		BGP+3707_2
	57 mm		BGP+5707_2
	77 mm		BGP+7707_2
8 mm	27 mm	8 F	BGP+2708_2
	37 mm		BGP+3708_2
	57 mm		BGP+5708_2
	77 mm		BGP+7708_2
9 mm	27 mm	8/9/10 F ⁵	BGP+2709_2
	37 mm		BGP+3709_2
	57 mm		BGP+5709_2
	77 mm		BGP+7709_2
10 mm	27 mm	8/9/10 F ⁵	BGP+2710_2
	37 mm		BGP+3710_2
	57 mm		BGP+5710_2
	77 mm		BGP+7710_2

For sizes with BEVAR indication, please consult the Instructions for Use.

All information provided in this brochure is not intended to be professional medical advice. The Instructions for Use (IFU) for the product must therefore be consulted before use of any product. Please be aware that the product and/or an indication might not be approved in your country, e.g. US. Please contact your Bentley representative for further information.

All data on file at Bentley InnoMed GmbH.

Be innovative

We focus on the development of unique products and strive to find creative answers to questions posed by the healthcare industry

Be dedicated

A marked customer focus and compliance to international standards help us to uphold unsurpassed quality in products and services

Be responsible

We care for better health, our region and the environment. Every member of our team is motivated to make a vital contribution to the quality of life

Be effective

We strive for time- and cost efficient solutions to satisfy the versatile interests of our customers and patients worldwide



Bentley InnoMed GmbH
 Lotzenäcker 3, 72379 Hechingen, Germany
 tel +49 7471 984 995 10
 mail info@bentley.global

CE 0123

© 2025 Bentley InnoMed GmbH | Germany
 All Rights Reserved | Art.-Nr. 408126-001/002