

PRODUCT DATA SHEET

DUAL BULB PASTETTE®


DUAL-BULB PASTETTES ARE DESIGNED TO DISPENSE AN EXACT VOLUME OF LIQUID WITH FOOLPROOF SIMPLICITY.


A quick and easy alternative to pipettes and tips they are particularly useful in self-test and Point of Care (POC) diagnostic kits where users may not have experience of exact volume pipetting.

Simply squeeze the upper bulb to fill, and squeeze again to dispense the exact volume. These pipettes are available unsterilised in volumes from 20-300 μ l.

- Fast measured transfers
- Easy to use
- Innovative
- Foolproof simplicity for exact volume delivery

**samples
available**

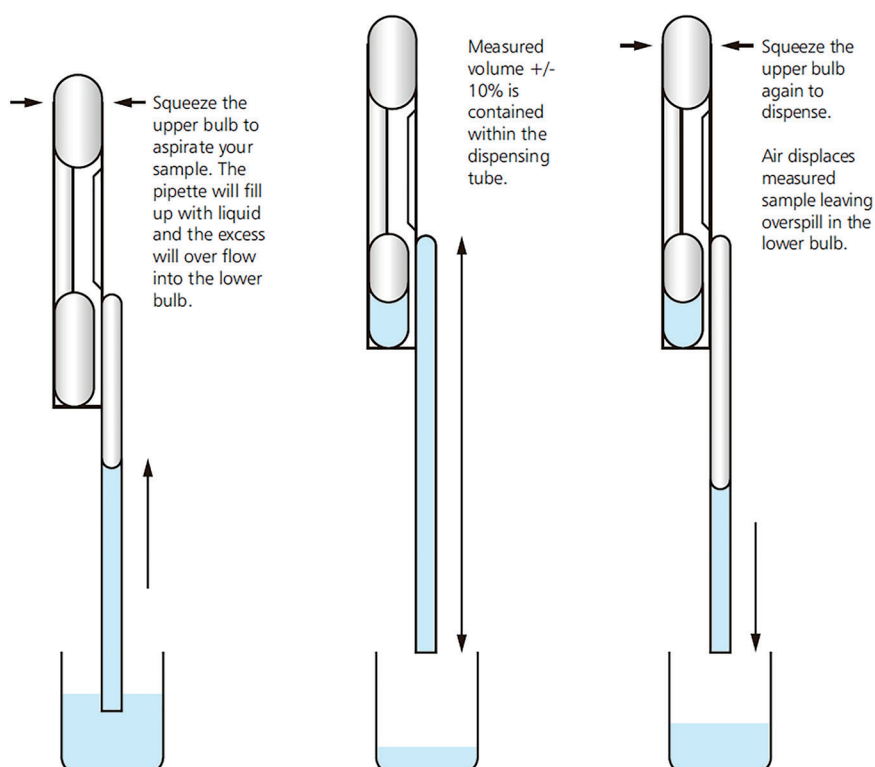

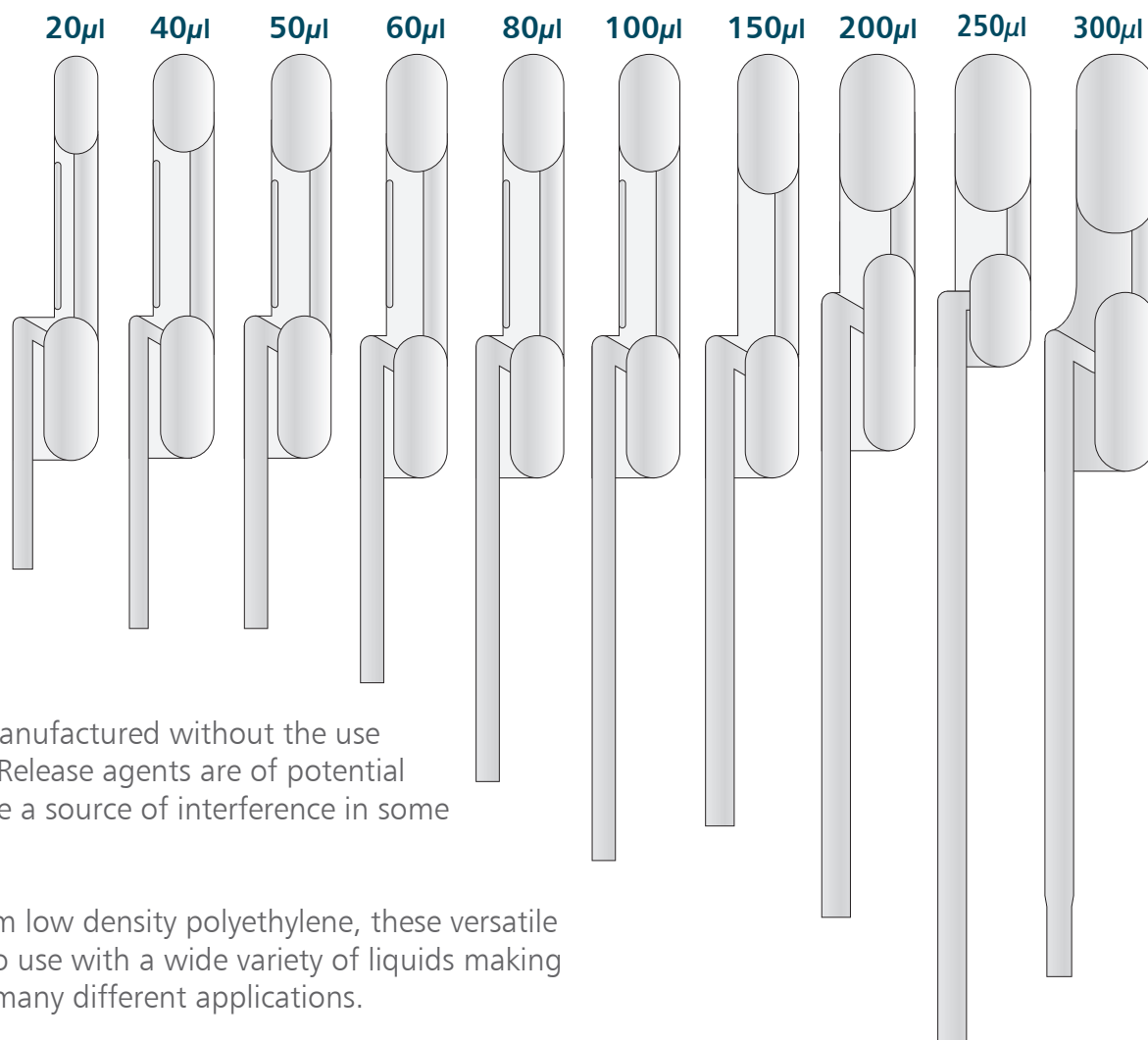

All Pastettes are manufactured without the use of release agents. Release agents are of potential concern as they are a source of interference in some biological assays.

Manufactured from low density polyethylene, these versatile pipettes are safe to use with a wide variety of liquids making them suitable for many different applications.

Dual-Bulb Pastette® (Non-Sterile)

| Cat. No. | Description | Bulb Draw | Length | Dispense Volume | Pack Size |
|----------|---------------------|-----------|--------|-----------------|-----------|
| LW4730 | Dual-Bulb Pastette® | 60µl | 53mm | 20µl | 5,000 |
| LW4740 | Dual-Bulb Pastette® | 110µl | 59mm | 40µl | 5,000 |
| LW4745 | Dual-Bulb Pastette® | 110µl | 59mm | 50µl | 5,000 |
| LW4750 | Dual-Bulb Pastette® | 146µl | 65mm | 60µl | 5,000 |
| LW4760 | Dual-Bulb Pastette® | 146µl | 75mm | 70µl | 5,000 |
| LW4770 | Dual-Bulb Pastette® | 146µl | 83mm | 100µl | 5,000 |
| LW4780 | Dual-Bulb Pastette® | 180µl | 80mm | 150µl | 5,000 |
| LW4785 | Dual-Bulb Pastette® | 275µl | 89mm | 200µl | 5,000 |
| LW4790 | Dual-Bulb Pastette® | 324µl | 102mm | 250µl | 5,000 |
| LW4795 | Dual-Bulb Pastette® | 385µl | 96mm | 300µl | 5,000 |

Certain volumes are also available in smaller pack sizes

| | | | | | |
|------------|---------------------|-------|-------|-------|-----|
| LW4730-500 | Dual-Bulb Pastette® | 60µl | 53mm | 20µl | 500 |
| LW4745-500 | Dual-Bulb Pastette® | 110µl | 59mm | 50µl | 500 |
| LW4770-500 | Dual-Bulb Pastette® | 146µl | 83mm | 100µl | 500 |
| LW4780-500 | Dual-Bulb Pastette® | 180µl | 80mm | 150µl | 500 |
| LW4790-500 | Dual-Bulb Pastette® | 324µl | 102mm | 250µl | 500 |